Assay Performance Characteristics:

Standard range: 40-0.08ng/mL Limit of Detection: 0.16ng/mL Background: OD<0.08 at 450nm

Coefficient of Determination: R-squared>0.98

Plate Template:

	1	2	3	4	5	6	7	8	9	10	11	12
A												
В												
С												
D												
Е												
F												
G												
Н												

References:











Indoor Biotechnologies, Inc. 700 Harris Street, Charlottesville, VA 22903 **United States**

Indoor Biotechnologies Ltd. Vision Court Caxton Place Cardiff. Wales CF23 8HA United Kingdom

Tel: (434) 984-2304 Fax:(434) 984-2709 Email: mail@inbio.com

Tel: +44 (0) 29 2167 4640 Email: Info@indoorbiotech.co.uk Indoor Biotechnologies, India Private Limited Bangalore Bioinnovation Centre, BioTech Park, Electronic City Phase 1 Bangalore-560100, India


Tel: +91-9901722009 Email: info@inbioindia.com



Ana o 3 ELISA 2.0 **Pre-coated Plate Kit**

Product Code: EPC-AO3-X

Sample curve:

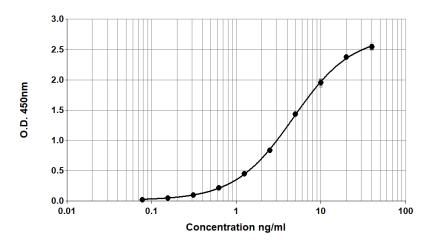

Contents:

Microtiter plate coated with anti-Ana o 3 monoclonal antibody 1H4

Ana o 3 allergen standard (white cap) Concentration: 400ng/mL

Biotinylated monoclonal antibody 4A11 (brown cap)

Streptavidin-peroxidase (blue cap)

Wash buffer (10x concentrate) Assay buffer (10x concentrate) TMB developing substrate Stop solution (0.5N sulfuric acid)

Store kit at 2-8°C Expiry: xxxxx xx, xxxx

> For research and commercial use in vitro: not for human in vivo or therapeutic use. An InBio[™] product



Certificate of Analysis

Pre-coated Plate: 96-well polystyrene microtiter plate coated with

monoclonal antibody 1H4 and treated with stabilizing

agent. Sealed in foil pouch with desiccant.

Monoclonal Antibody: 1H4

Immunogen: Ana o 3

Isotype:

Specificity: Binds to a specific epitope on Anacardium occidentale

allergen Ana o 3

Purification: Produced in cell culture and purified by affinity

chromatography using Protein A. Single heavy and light

chain bands on SDS-PAGE.

Lot Number: xxxxx

Detection Antibody: 4A11

Immunogen: Ana o 3

Isotype: Mouse IgG2b/kappa

Specificity: Binds to a specific epitope on *Anacardium occidentale*

allergen Ana o 3

Purification: Produced in cell culture and purified by affinity

chromatography using Protein A. Single heavy and

light chain bands on SDS-PAGE.

Biotinylation: Biotinylated and titrated for use in ELISA at 1/1000

dilution. Prepared in 1% BSA/50% glycerol/PBS, pH 7.4,

0.22µm filtered, preservative free.

Lot Number: xxxxx

Allergen Standard: Purified recombinant Ana o 3 prepared in 1% BSA/50%

glycerol/PBS, pH 7.4.

Concentration: 400ng/mL (based on amino acid analysis)

Lot Number: xxxxx

Materials required, but not provided:

- Type I ultrapure water or 18.2MΩ de-ionized water
- Volumetric measuring equipment (e.g. serological pipettes, graduated cylinders)
- Clean containers for buffer and reagent preparation
- Calibrated single and multi-channel micropipettes and tips
- Vortex mixer
- Plate reader capable of reading absorbance at 450nm
- Analysis software (recommended, but not required)

Protocol

Please read the entire protocol before starting the assay

Bring all reagents to room temperature before use

1. Prepare 1x working dilutions of the 10x wash and assay buffers in clean containers using $18.2M\Omega$ de-ionized water or Type I ultrapure water. For one plate:

Wash buffer: add 15mL concentrate to 135mL water Assay buffer: add 2.5mL concentrate to 22.5mL water

Adjust volumes accordingly for multi-plate assays. Diluted buffers may be stored at 4 C for up to 1 week

*Use of an automated plate washer will require a larger volume of wash buffer than provided. Most PBS-based wash buffers should be compatible with the assay. Visit the support page at www.inbio.com for recommended wash buffer formulary.

- 2. Remove the plate from the foil pouch and wash by adding 150µL wash buffer to each well. Empty the wells by inverting the plate and then tap on absorbent paper to remove residual buffer. Repeat the wash cycle two more times.
- 3. Add standards, samples, and blanks to the plate (final volume in all wells is 100µL).

Standards: add 180 μ L assay buffer into wells A1 and B1, and 100 μ L into remaining wells of rows A and B. Vortex the Ana o 3 standard and add 20 μ L to wells A1 and B1. Mix well by pipetting up and down 7-10 times and then transfer 100 μ L into wells A2 and B2. Mix well and continue the serial doubling dilution scheme across the plate to column 10.

The assay buffer in wells A11, B11 and A12, B12 will serve as **Blanks**.

Samples: dust extracts are routinely tested starting at 1/10 dilution and can be prepared directly on the pre-coated plate: add 20µL sample to 180µL assay buffer. Mix, then transfer 100µL into 100µL assay buffer in the next well. Continue across the plate for the desired number of dilutions. A minimum of three dilutions per sample should be tested; 6-12 dilutions are recommended.

Air filter extracts, allergen extracts, and other types of samples may require a different dilution scheme.

*Sample dilutions may also be prepared in tubes or on a 96-well dilution plate and transferred to the pre-coated plate.

- 4. Incubate the plate at room temperature (away from direct sunlight) for 1 hour.
- Wash the plate 3x with 150µL wash buffer per well. Vortex the biotinylated 4A11 and prepare a 1:1,000 detection antibody/conjugate mix by adding 10µL biotinylated 4A11and 10µL streptavidin-peroxidase to 10mL assay buffer.
 Mix thoroughly and add 100µL to each well.
- 6. Incubate the plate at room temperature (away from direct sunlight) for 1 hour.
- 7. Pour the TMB substrate and stop solution into separate basins so they are ready to use in the next step. Wash the plate 3x with 150µL wash buffer per well.
- Use a <u>multi-channel</u> pipette to add 100μL TMB to each well. Gently tap the plate and monitor the reaction as the blue color develops. Once OD450 reaches 0.08-0.09 for Standard 1, use a <u>multi-channel</u> pipette to add 50μL stop solution to each well (the color will change to yellow).
- 9. Read the plate at 450nm. The OD for Standard 1 should be between 1.2 and 3.5, with an ideal range of 2.0 2.5.

A list of frequently asked questions and troubleshooting guide can be found under the 'Support' tab on our web site: www.inbio.com.