Assay Performance Characteristics:

Standard range: 100-0.02ng/mL Limit of Detection: 0.78ng/mL Background: OD<0.08 at 450nm Coefficient of Determination: R-squared>0.98

Plate Template:

| | 1 | 2 | 3 | 4 | 5 | 6 | 7 | 8 | 9 | 10 | 11 | 12 |
|---|---|---|---|---|---|---|---|---|---|----|----|----|
| А | | | | | | | | | | | | |
| В | | | | | | | | | | | | |
| С | | | | | | | | | | | | |
| D | | | | | | | | | | | | |
| Е | | | | | | | | | | | | |
| F | | | | | | | | | | | | |
| G | | | | | | | | | | | | |
| Н | | | | | | | | | | | | |

References:

- Luczynska CM, Arruda LK, Platts-Mills TA, Miller JD, Lopez M, Chapman MD. A two-site monoclonal antibody ELISA for the quantification of the major Dermatophagoides spp. allergens, Der p I and Der f I. J Immunol Methods 1989; 118(2):227-235.
- Custovic A, Taggart SC, Francis HC, Chapman MD, Woodcock A. Exposure to house dust mite allergens and the clinical activity of asthma. J Allergy Clin Immunol 1996; 98 (1):64-72.



Indoor Biotechnologies, Inc. 1216 Harris Street, Charlottesville,

Charlottesville, VA 22903 United States

Tel: (434) 984-2304 Fax:(434) 984-2709 E-mail: mail@inbio.com



*i*es[™] **●**

Indoor Biotechnologies Ltd. Cardiff Medicentre Heath Park Cardiff CF14 4UJ Wales

Tel: +44 (0) 2920 68 2115 Fax: +44 (0) 2920 68 2044 E-mail: info@indoorbiotech.co.uk


RG



Der p 1 ELISA 2.0 Pre-coated Plate Kit

Product Code: EPC-DP1-X



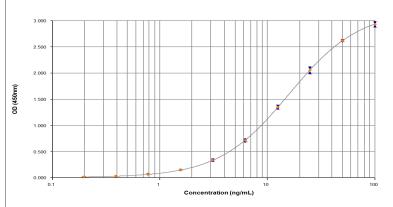

Contents:

Plate: Pre-coated with anti-Der p 1 monoclonal antibody 5H8

- Vial 1: (white top) Der p 1 allergen standard Concentration: 1000ng/ml
- Vial 2: (brown) Biotinylated monoclonal antibody 4C1
- Vial 3: (red top) Streptavidin-peroxidase

Bottle 1: Wash buffer, (10x concentrate) Bottle 2: Assay buffer, (10x concentrate) Bottle 3: TMB developing substrate Bottle 4: Stop solution (0.5N sulfuric acid)

Store kit at 2-8^oC Expiry: 6 months from date of receipt

> For research and commercial use in vitro: not for human in vivo or therapeutic use.

INDOOR and ne Registered Trademarks of INDOOR Biotechnologies, Ltd., Manchester U.K.

| | Certificate of Analysis | | | | |
|--|--|--|--|--|--|
| Pre-coated Plate: | 96-well polystyrene microtiter plate coated with monoclonal antibody 5H8 and treated with stabilizing agent. Sealed in foil pouch with desiccant. | | | | |
| Monoclonal Antibody: Immunogen: Isotype: Specificity: Purification: Lot Number: | 5H8 (clone 5H8 C12 D8) Der p 1 Mouse IgG2A Binds to an epitope on dust mite <i>Dermatophagoides</i> <i>pteronyssinus</i> allergen, Der p 1. Produced in tissue culture and purified by affinity chromatography using Protein A. Single heavy and light chain bands on SDS-PAGE. XXXXX | | | | |
| Detection Antibody: | 4C1 (clone 4C1 B8 3F8) | | | | |
| Immunogen: Isotype: Specificity: Purification: Biotinylation: Lot Number: | Der f 1 Mouse IgG1 Binds to an epitope on dust mite <i>Dermatophagoides</i> Group 1 allergens (Der f 1, Der p 1, Der m 1, Eur m 1). Produced in ascites and purified by affinity chromatography using Protein A. Single heavy and light chain bands on SDS-PAGE. Biotinylated and titrated for use in ELISA at 1/1000 dilution. Prepared in 1% BSA/50% glycerol/PBS, pH 7.4, 0.22µm filtered, preservative free. XXXXX | | | | |
| Allergen Standard: | Purified natural Der p 1 prepared in 1% BSA/50% glycerol/PBS, pH 7.4. | | | | |
| Concentration: Lot Number: | 1,000ng/mL (based on amino acid analysis) XXXXX | | | | |

Materials required, but not provided:

- Type I ultrapure water or $18.2 M\Omega$ de-ionized water
- Volumetric measuring equipment (e.g. serological pipette, graduated cylinder)
- Clean containers for buffer and reagent preparation
- Calibrated single and multi-channel micropipettes and tips
- Vortex mixer
- Plate reader capable of reading absorbance at 450nm
- Analysis software (recommended, but not required)

Protocol

Please read entire protocol before starting the assay

Bring all reagents to room temperature and vortex before use.

- 1. Allow the pre-coated plate to reach room temperature while in the sealed pouch.
- Prepare 1x working dilutions of the 10x wash and assay buffers in clean containers using 18.2MΩ de-ionized water or Type I ultrapure water. For one plate:
 Wash buffer: add 15mL concentrate to 135mL water

Assay buffer: add 2.5mL concentrate to 22.5mL water Adjust volumes accordingly for multi-plate assays. *Diluted buffers may be stored at 4^oC for up to 1 week

- Remove the plate from the foil pouch and wash by adding 150µL wash buffer to each well. Empty the wells by inverting the plate and tapping on absorbent paper. Repeat the wash 2x.
- Add standards, samples and blanks to the plate (final volume in all wells should be 100µL).
 Standards: pipette 180µL assay buffer into wells A1 and B1 and 100µL into remaining wells of rows A and B. Add 20µL Der p 1 standard to wells A1 and B1. Mix well and transfer 100µl into wells A2 and B2. Continue across the plate to wells A10 and B10 to make 10 serial doubling dilutions.
 Samples: dust extracts for Der p 1 analysis are routinely diluted two-fold starting at 1/10. Other types of samples, like air filter extracts and allergen extracts, may require different dilutions. It is recommended to test each sample at a minimum of three dilutions; 6-12 are recommended.
 Blanks: add assay buffer to wells A11, B11 and A12, B12.
- 5. Return the plate to the foil pouch or cover with plate sealer and incubate for 1 hour at room temperature.
- Wash wells 3x with 150µL wash buffer. Prepare a 1:1,000 detection antibody/conjugate mix by adding 10µL biotinylated 4C1 and 10µL streptavidin-peroxidase to 10mL assay buffer. Mix well and add 100µl to each well. Cover and incubate for 1 hour at room temperature.
- 7. Wash wells 3x with 150µL wash buffer. Add 100µl TMB substrate to each well and monitor the reaction as the blue color develops. When OD450nm reaches 0.08 for the first standard in wells A1 and B1 (generally within 1-5 minutes), add 50µL stop solution (the color will change to yellow).
- 8. Read the plate at 450nm. The ideal OD for standard 1 is 2.0-2.5.

Notes:

The allergen standard is recommended for immunoassay calibration purposes only.

A list of frequently asked questions and troubleshooting guide can be found under the 'Support' tab on our web site: www.inbio.com.