


Equ c 4 ELISA kit (103/14G4)

Product Code: EL-EC4

Lot Number: xxxxx

Sample Curve:

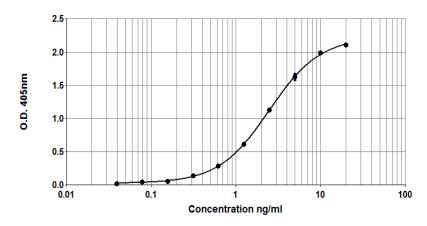

Content:

- Vial 1 (red top) 100 µL Monoclonal antibody 103 Concentration: 1mg/ml in PBS
- Vial 2 (white top) 400 µL Equ c 4 Standard Concentration: 200U/ml Equ c 4
- Vial 3 (brown) 100 µL Biotinylated monoclonal antibody 14G4 Dilute: 1:1000 for use

Storage: The ELISA kit should be stored at 4°C

For research and commercial use in vitro: not for human in vivo or therapeutic use.



www.inbio.com

Indoor Biotechnologies, Inc. 1216 Harris Street, Charlottesville, VA 22903 United States

Tel: (434) 984-2304 Fax:(434) 984-2709 E-mail: mail@inbio.com Indoor Biotechnologies Ltd. Cardiff Medicentre Heath Park Cardiff CF14 4UJ Wales, United Kingdom

Tel: +44 (0) 2920 68 2115 Fax: +44 (0) 2920 68 2044 E-mail: info@indoorbiotech.co.uk

Certificate of Analysis

Monoclonal Antibody:	103
Immunogen:	Horse Allergen
Specificity:	Binds to an epitope on horse Equus caballus
	allergen, Equ c 4.
Concentration:	1.0mg/ml in phosphate buffered saline. Based
	on A280 for IgG (1.42=1mg/ml) 0.22µm filtered,
	0.02% sodium azide.
Lot Number:	XXXXX

Monoclonal Antibody:	14G4
Immunogen:	Horse Allergen
Specificity:	Binds to an epitope on horse <i>Equus caballus</i> allergen, Equ c 4.
Concentration:	1.0mg/ml in phosphate buffered saline. Based on A280 for IgG (1.42=1mg/ml) 0.22µm filtered, 0.02% sodium azide.
Biotinylation:	Biotinylated and titrated for use in ELISA at 1/1000 dilution.
Lot Number:	XXXXX

Allergen Standard:Equ c 4 StandardComposition:Horse hair extract prepared in 1% BSA/50%
glycerol/PBS, pH 7.4Concentration:200U/mlLot Numberxxxxx

ELISA Protocol for Equ c 4.

- 1. Coat polystyrene microtiter plates (NUNC Maxisorp Cert. NUNC catalog # 439454) with 100µl mAb 103 at 10µl/10ml, i.e. 1/1000 dilution of stock, in 50mM carbonate-bicarbonate buffer, pH 9.6, incubate overnight at 4°C.
- Wash wells 3x with PBS-0.05% Tween 20, pH 7.4 (PBS-T). Incubate for 30 min. at room temperature with 100µl/well of 1% BSA, PBS-T. Wash 3x with PBS-T.
- 3. Use doubling dilutions of the Horse Standard to make a control curve ranging from 20 0.04U/ml Equ c 4: Pipette 20µl standard into 180µl 1% BSA, PBS-T into wells A1 and B1 on the ELISA plate. Mix well and transfer 100µl across the plate into 100µl 1% BSA, PBS-T diluent to make 10 serial doubling dilutions. Wells A11, B11 and A12, B12 should contain only 1% BSA, PBS-T as blanks.
- 4. Add 100µl of diluted allergen samples and incubate for 1 hour at room temperature. House dust extracts for Equ c 4 analysis are routinely diluted two-fold from1/10-1/80. Other sample types, like air filter extracts and allergen extracts, may require different dilutions.
- Wash wells 3x with PBS-T and add 100µl diluted biotinylated anti-Equ c 4 mAb 14G4. The antibody solution should be diluted 1/1000 in 1% BSA, PBS-T. Incubate for 1 hour at room temperature.
- Wash wells 3x with PBS-T and add 100μl diluted Streptavidin-Peroxidase (Sigma S5512, 0.25mg reconstituted in 1ml distilled water). The reconstituted Streptavidin should be diluted 1/1000 in 1% BSA, PBS-T. Incubate for 30 minutes at room temperature.
- 7. Wash wells 3x and develop the assays by adding 100μ I 1mM ABTS in 70mM citrate phosphate buffer, pH 4.2 and 1/1000 dilution of H₂O₂. Read the plate when the absorbance at 405nm reaches 2.0-2.4.

Notes:

The Horse Allergen Standard is recommended for immunoassay calibration purposes only. Not recommended for in-vitro antibody measurements, T cell studies, immunization purposes, or other uses.

Buffer recipes, storage conditions and a list of frequently asked questions can be found under "Protocols" on our web site: www.inbio.com.

For research and commercial use in vitro: not for human in vivo or therapeutic use.