

Product Datasheet

VEGF R2/KDR/FIk-1 Antibody

NBP2-36428SS

Unit Size: 0.025 mg

Store at 4C short term. Aliquot and store at -20C long term. Avoid freeze-thaw cycles.

www.novusbio.com



support@novusbio.com

Protocols, Publications, Related Products, Reviews, Research Tools and Images at:
www.novusbio.com/NBP2-36428

Updated 9/28/2014 v.20.1

NBP2-36428SS

VEGF R2/KDR/FIk-1 Antibody (1B6)


Product Information	
Unit Size	0.025 mg
Concentration	1.0 mg/ml
Storage	Store at 4C short term. Aliquot and store at -20C long term. Avoid freeze-thaw cycles.
Clonality	Monoclonal
Clone	1B6
Preservative	0.05% Sodium Azide
Isotype	IgG1 Kappa
Purity	Protein G purified
Buffer	PBS

Product Description	
Host	Mouse
Gene ID	3791
Gene Symbol	KDR
Species	Human
Immunogen	Partial recombinant human VEGFR2 protein (between amino acids 300-600) [UniProt P35968]


Product Application Details	
Applications	Western Blot, Flow (Cell Surface), Immunocytochemistry/Immunofluorescence, Immunohistochemistry, Immunohistochemistry-Paraffin
Recommended Dilutions	Flow (Cell Surface), Immunocytochemistry/Immunofluorescence 1:50 - 1:100, Immunohistochemistry, Immunohistochemistry-Paraffin 5 ug/ml, Western Blot 1 ug/ml

Images


Western Blot: VEGF Receptor 2 Antibody (1B6) [NBP2-36428] - Western blot analysis of VEGF Receptor 2 (1B6) in purified human recombinant protein

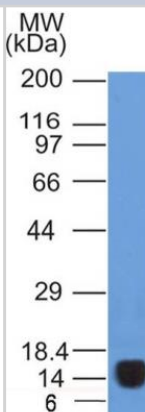


Immunocytochemistry/Immunofluorescence: VEGF Receptor 2 Antibody (1B6) [NBP2-36428] - ICC detection of VEGF Receptor 2 in HUVEC cells at 1:50

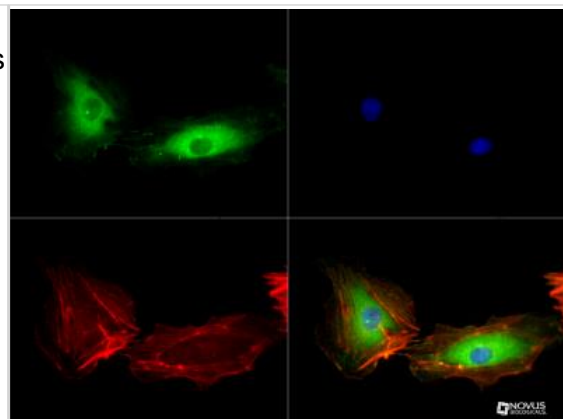


Immunohistochemistry-Paraffin: VEGF Receptor 2 Antibody (1B6) [NBP2-36428] - IHC analysis of formalin-fixed paraffin-embedded tissue section of a malignant stromal tumor of the human small bowel with mouse monoclonal VEGF Receptor 2 antibody (clone 1B6) at 5 ug/ml concentration. The cancer cells depicted a specific and strong staining in the cytoplasm.

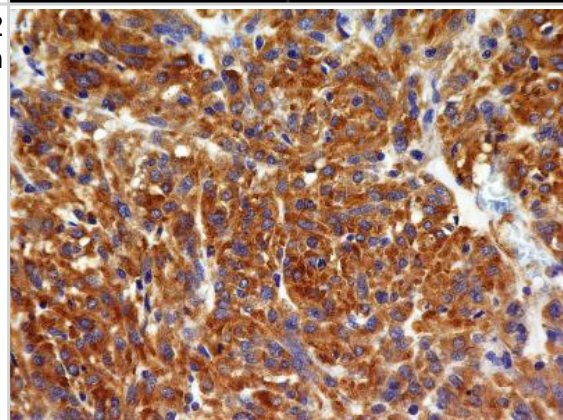


Flow Cytometry: VEGF Receptor 2 Antibody (1B6) [NBP2-36428] - HepG2 cells were stained (surface) with VEGFR2 antibody (clone: 2B5.1B6; orange) or isotype control (mouse IgG1; blue) and positive staining observed using FITC conjugated mouse anti-IgG(H+L) secondary antibody. Live cells (PPI negative) were gated for analysis

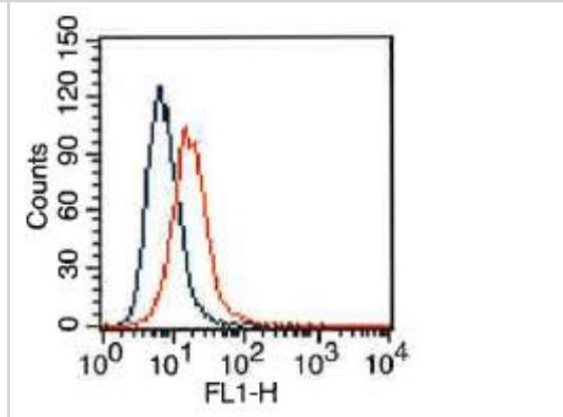


Immunohistochemistry-Paraffin: VEGF Receptor 2 Antibody (1B6) [NBP2-36428] - IHC analysis of formalin-fixed paraffin-embedded tissue section of human normal kidney with mouse monoclonal VEGF Receptor 2 antibody (clone 1B6) at 5 ug/ml concentration. The cuboidal epithelial cells of various tubules/ducts showed a specific cytoplasmic immunopositivity of VEGFR2 protein.

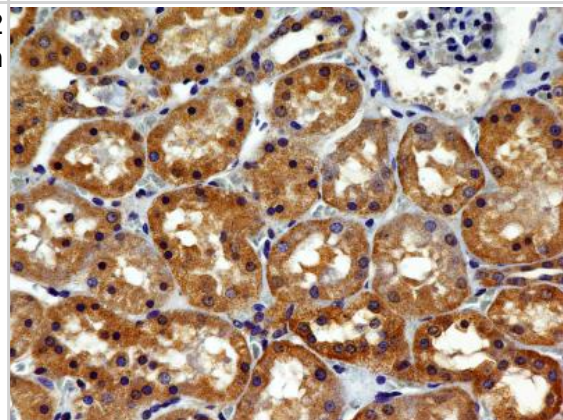

Immunohistochemistry-Paraffin: VEGF Receptor 2 Antibody (1B6) [NBP2-36428] - IHC analysis of formalin-fixed paraffin-embedded transverse section of normal human skeletal muscle with mouse monoclonal VEGF Receptor 2 antibody (clone 1B6) at 5 ug/ml concentration. The myocytes showed a moderate diffused to granular staining in the cytoplasm.

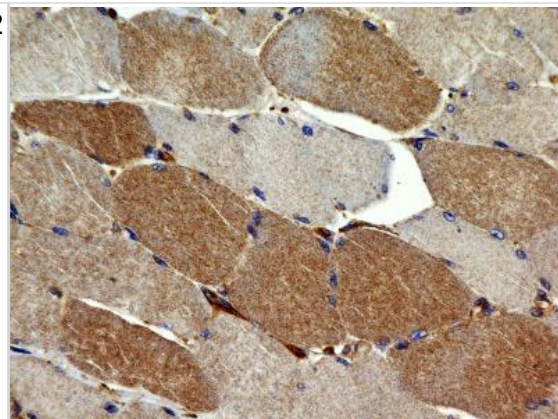

Immunohistochemistry-Paraffin: VEGF Receptor 2 Antibody (1B6) [NBP2-36428] - IHC analysis of formalin-fixed paraffin-embedded tissue section of human liver cancer with mouse monoclonal VEGF Receptor 2 antibody (clone 1B6) at 5 ug/ml concentration. The cancer cells showed a specific and strong staining in the cytoplasm.

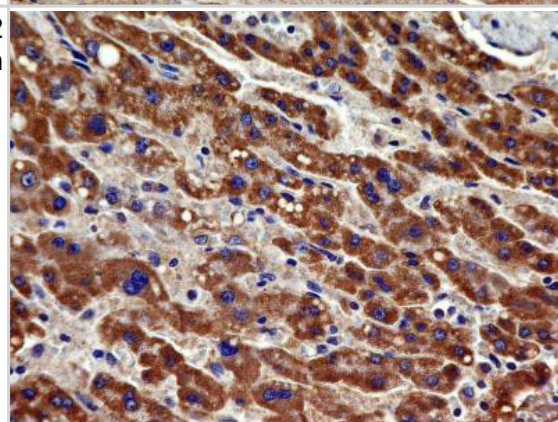

Procedures

Western Blot Protocol for VEGF Receptor 2 Antibody (NBP2-36428)

Western Blot Protocol

1. Perform SDS-PAGE on samples to be analyzed, loading 25 ug of total protein per lane.
2. Transfer proteins to membrane according to the instructions provided by the manufacturer of the membrane and transfer apparatus.
3. Stain according to standard Ponceau S procedure (or similar product) to assess transfer success, and mark molecular weight standards where appropriate.
4. Rinse the blot.
5. Block the membrane using standard blocking buffer for at least 1 hour.
6. Wash the membrane in wash buffer three times for 10 minutes each.
7. Dilute anti-HDAC5 primary antibody in blocking buffer and incubate 1 hour at room temperature.
8. Wash the membrane in wash buffer three times for 10 minutes each.
9. Apply the diluted HRP conjugated secondary antibody in blocking buffer (as per manufacturers instructions) and incubate 1 hour at room temperature.
10. Wash the blot in wash buffer three times for 10 minutes each (this step can be repeated as required to reduce background).
11. Apply the detection reagent of choice in accordance with the manufacturers instructions.

Note: Tween-20 can be added to the blocking or antibody dilution buffer at a final concentration of 0.05-0.2%.



Immunohistochemistry-Paraffin Protocol for VEGF Receptor 2 Antibody (NBP2-36428)

1. Deparaffinize the tissue sections by immersing the slides in Xylene with two changes for 10 min each. Sections should not get dried at any stage from this point.
2. Rehydrate the tissue sections by immersing the slides in decreasing grades of ethanol as follows:
 - a. Immerse in 100% ethanol with 2 changes for 5 minutes each
 - b. Immerse in 95% ethanol with 2 changes for 5 minutes each
 - c. Immerse in 90% ethanol for 5 minutes
 - d. Immerse in 70% ethanol for 5 minutes
 - e. Immerse in 50% ethanol for 5 minutes
 - f. Immerse in distilled water for 5 minutes
3. Antigen Retrieval (Microwave Method):
 - a. Immerse the slides in a microwave compatible tray containing 10 mM Sodium Citrate buffer (pH 6.0) with 0.05% Tween 20.
 - b. Boil the slides and maintain the sub-boiling temperature for 5 minutes in the microwave. Thereafter, take out the tray very carefully and cool it at room temperature (RT) for about 30 minutes.
 - c. Wash the slides 3 times, 3 minutes each by immersing them in TBST (Tris Buffered Saline having 0.05% Tween 20).
4. Quenching of Endogenous Peroxidase:
 - a. Incubate the slides in 3% hydrogen peroxide prepared in methanol for 15 minutes (at RT, in dark conditions).
 - b. Wash the slides in TBST 3 times, 3 minutes each.
5. Protein Blocking:
 - a. Incubate the sections with background sniper solution at RT for 15 minutes (Biocare Medicals, USA).
 - b. Wash the sections 3 times, 3 min each by immersing the slides in TBST.
6. Primary Antibody:
 - a. Dilute the primary antibody at 5ug/ml concentration using PBS as a diluent.
 - b. Incubate the sections with diluted primary antibody for 90 minutes at RT in a humidified chamber.
 - c. Thereafter, wash the slides 4 times, 5 minutes each with TBST.
7. Probe (Secondary Reagent):
 - a. Incubate with MACH 1 Mouse probe for 15 minutes at RT.
 - b. Incubate for 30 min at room temperature with HRP-Polymer (Biocare Medical, USA).
 - c. Wash the slides with TBST 4 times, 5 minutes each
8. Chromogen:
 - a. Mix 32ul of DAB Chromogen with 1 ml of DAB substrate buffer (Biocare Medical, USA).
 - b. Apply 200ul DAB mixture/section and incubate at RT in dark conditions (few seconds - 5 minutes).
 - c. As soon as an appropriate color develops, rinse the slides with deionized water (2-3 brief rinses).
9. Counter stain:
 - a. Counter stain with Hematoxylin for 30 seconds (Vector Labs, USA).
 - b. Wash in deionized water for 1-2 minutes to clear the extra stain.
 - c. Incubate the slides in bluing solution or Scott's water twice for 2 minutes each time.
10. Dehydrate the sections in increasing grades of alcohols:
 - a. 50% alcohol for 1 minute
 - b. 70% for 1 minute
 - c. 90% for 1 minute
 - d. 95% for 1 minute
 - e. 100% for 1 minute
 - f. Xylene with 2 changes for 2 minutes each
11. Mount with DPX mount and cover-slip glass (Fisher Scientific, USA), carefully not allowing any air bubbles to enter.

NOTE:- This protocol is provided as a reference tool only. Depending upon the type of tissues /tissue processing and reagents employed, the end user will need to optimize the final conditions for achieving an expected staining.



Immunocytochemistry/Immunofluorescence Protocol for VEGF Receptor 2 Antibody (NBP2-36428)

Immunocytochemistry Protocol

Culture cells to appropriate density in 35 mm culture dishes or 6-well plates.

1. Remove culture medium and add 10% formalin to the dish. Fix at room temperature for 30 minutes.
2. Remove the formalin and add ice cold methanol. Incubate for 5-10 minutes.
3. Remove methanol and add washing solution (i.e. PBS). Be sure to not let the specimen dry out. Wash three times for 10 minutes.
4. To block nonspecific antibody binding incubate in 10% normal goat serum from 1 hour to overnight at room temperature.
5. Add primary antibody at appropriate dilution and incubate at room temperature from 2 hours to overnight at room temperature.
6. Remove primary antibody and replace with washing solution. Wash three times for 10 minutes.
7. Add secondary antibody at appropriate dilution. Incubate for 1 hour at room temperature.
8. Remove antibody and replace with wash solution, then wash for 10 minutes. Add Hoechst 33258 to wash solution at 1:25,000 and incubate for 10 minutes. Wash a third time for 10 minutes.
9. Cells can be viewed directly after washing. The plates can also be stored in PBS containing Azide covered in Parafilm (TM). Cells can also be cover-slipped using Fluoromount, with appropriate sealing.

*The above information is only intended as a guide. The researcher should determine what protocol best meets their needs. Please follow safe laboratory procedures.





Novus Biologicals USA

8100 Southpark Way, A-8
Littleton, CO 80120
USA
Phone: 303.730.1950
Toll Free: 1.888.506.6887
Fax: 303.730.1966
novus@novusbio.com

Novus Biologicals Canada

461 North Service Road West, Unit B37
Oakville, ON L6M 2V5
Canada
Phone: 905.827.6400
Toll Free: 855.668.8722
Fax: 905.827.6402
canada@novusbio.com

Novus Biologicals Europe

12 Cambridge Science Park
Cambridge, CB4 0FQ
United Kingdom
Phone: +44 (0)1223 426001
Fax: +44 (0)871 971 1635
europe@novusbio.com

General Contact Information

www.novusbio.com
Technical Support: technical@novusbio.com
Orders: orders@novusbio.com
General: novus@novusbio.com

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Primary Antibodies are guaranteed for 1 year from date of receipt.

For more information on our guarantee, please visit www.novusbio.com/guarantee.

