

ELISpot

Mouse TNF- α

Catalog Number EL410

For the quantitative determination of the frequency of cells releasing mouse TNF- α .

This package insert must be read in its entirety before using this product.

**FOR RESEARCH USE ONLY.
NOT FOR USE IN DIAGNOSTIC PROCEDURES.**

TABLE OF CONTENTS

Contents	Page
INTRODUCTION	2
PRINCIPLE OF THE ASSAY	3
ELISpot SCHEMATIC	4
LIMITATIONS OF THE PROCEDURE	5
PRECAUTIONS	5
MATERIALS PROVIDED.	5
STORAGE	5
OTHER SUPPLIES REQUIRED	6
TECHNICAL HINTS	6
REAGENT PREPARATION	6
SAMPLE PREPARATION	7
ASSAY PROCEDURE	7
CALCULATION OF RESULTS	8
REPRODUCIBILITY DATA.	8
TROUBLESHOOTING GUIDE	9
REFERENCES	10
ASSAY RECORD TEMPLATE	11

MANUFACTURED AND DISTRIBUTED BY:

R&D Systems, Inc.
614 McKinley Place NE
Minneapolis, MN 55413
United States of America

TELEPHONE: (800) 343-7475
(612) 379-2956
FAX: (612) 656-4400
E-MAIL: info@RnDSystems.com

DISTRIBUTED BY:

R&D Systems Europe, Ltd.
19 Barton Lane
Abingdon Science Park
Abingdon, OX14 3NB
United Kingdom

TELEPHONE: +44 (0)1235 529449
FAX: +44 (0)1235 533420
E-MAIL: info@RnDSystems.co.uk

R&D Systems China Co. Ltd.
24A1 Hua Min Empire Plaza
726 West Yan An Road
Shanghai PRC 200050

TELEPHONE: +86 (21) 52380373
FAX: +86 (21) 52371001
E-MAIL: info@RnDSystemsChina.com.cn

INTRODUCTION

Mouse TNF- α is a pleiotropic cytokine that plays a central role in inflammation and apoptosis (1 - 5). It is synthesized initially as a 26 kDa type II transmembrane (TM) glycoprotein that is 235 amino acids (aa) in length (6 - 8). It contains a 35 aa cytoplasmic domain, a 21 aa TM segment, and a 179 aa extracellular region that is divided into a 23 aa linker region, and a 156 aa mature molecule (6). TNF- α is assembled intra-cellularly to form a transmembrane, noncovalently linked homotrimeric protein. The 156 aa soluble form of trimeric TNF is released from the C-terminus of the TM protein through the activity of TACE, a membrane-bound disintegrin metalloprotease (9). Mouse TNF- α shares 80% aa identity with human TNF- α , and 95% aa identity with rat TNF- α (6, 10, 11). Mouse and human TNF- α show significant cross-species bioactivity (6). Mouse cells known to express TNF- α include neurons (12), CD8⁺ T cells (13), mast cells (14), macrophages (15), astrocytes, monocytes and microglia (16) and dendritic cells (17). Rat cells known to express TNF- α include endothelial cells (18) and pancreatic acinar cells (19).

There are two mouse receptors for TNF- α . TNF RI/TNFRSF1A is a 55 - 60 kDa type I TM glycoprotein that contains a death domain and binds TNF- α with a high affinity (20 - 22). TNF RII/TNFRSF1B is a 75 - 80 kDa TM glycoprotein that also binds TNF- α with high affinity (20). Although TNF RI and TNF RII share a general structural design, TNF RII has little amino acid identity to TNF RI. Unlike TNF RI, which is species indifferent, TNF RII is partially species specific and prefers mouse TNF- α (20).



TNF- α is involved in a number of pathophysiological processes. It is considered one of the prototypical pro-inflammatory molecules and is induced in macrophages by gram-negative bacteria (LPS). Once expressed, TNF- α is reported to promote inflammatory cell infiltration by upregulating leukocyte adhesion molecules on endothelial cells, serve as a chemotactic agent for monocytes, and activate phagocyte killing mechanisms (via increased NO₂⁻/O₂⁻/H₂O₂) (23). TNF- α also has documented antiviral activity, possibly acting through iNOS or IFN- β (24).

The Mouse TNF- α ELISpot assay is designed for the detection of mouse TNF- α secreting cells at the single cell level, and it can be used to quantitate the frequency of mouse TNF- α secreting cells. ELISpot assays are well suited for monitoring cytokine responses to various stimuli, treatments, and therapies, and have been used for the quantitation of various TNF- α secreting cell types. Other methods used for the assessment of TNF- α secreting cells, such as the chromium release assay with quantitation by limiting dilution, are tedious and require previous *in vitro* expansion of TNF- α producing cells for several days. These assays typically are not suitable for measuring infrequent TNF- α responses that occur at less than 1 in 1000. ELISpot assays are highly reproducible and sensitive and can be used to measure responses with frequencies well below 1 in 100,000. ELISpot assays do not require prior *in vitro* expansion of TNF- α producing cells, and thus are suitable for high-throughput analysis using only small volumes of primary cells. As such, ELISpot assays are useful tools for research into areas as diverse as inflammation, vaccine development, and the monitoring of various clinical trials.


PRINCIPLE OF THE ASSAY

The enzyme-linked immunospot (ELISpot) assay was originally developed for the detection of individual B cells secreting antigen-specific antibodies (25, 26). This method has since been adapted for the detection of individual cells secreting specific cytokines or other antigens (27, 28). ELISpot assays employ the quantitative sandwich enzyme-linked immunosorbent assay (ELISA) technique. A polyclonal antibody specific for mouse TNF- α has been pre-coated onto a PVDF (polyvinylidene difluoride)-backed microplate. Appropriately stimulated cells are pipetted into the wells and the microplate is placed into a humidified 37° C CO₂ incubator for a specified period of time. During this incubation period, the immobilized antibody in the immediate vicinity of the secreting cells binds secreted TNF- α . After washing away any cells and unbound substances, a biotinylated polyclonal antibody specific for mouse TNF- α is added to the wells. Following a wash to remove any unbound biotinylated antibody, alkaline-phosphatase conjugated to streptavidin is added. Unbound enzyme is subsequently removed by washing and a substrate solution (BCIP/NBT) is added. A blue-black colored precipitate forms and appears as spots at the sites of cytokine localization, with each individual spot representing an individual TNF- α secreting cell. The spots can be counted with an automated ELISpot reader system or manually using a stereomicroscope.


ELISpot SCHEMATIC

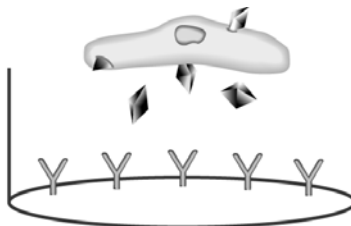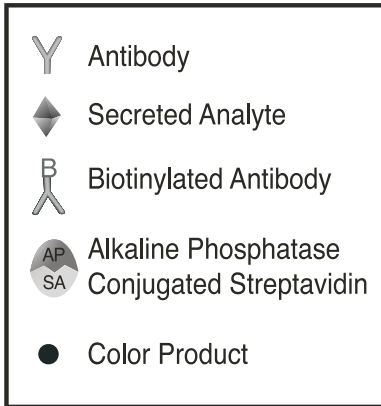


Incubate TNF- α -secreting cells in an antibody-coated well.

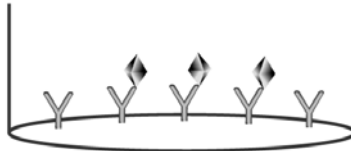


Remove cells by washing. Secreted TNF- α is captured by the immobilized antibody.

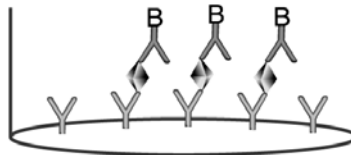


Incubate with biotinylated anti-TNF- α antibody.

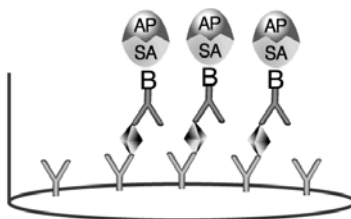

Incubate with alkaline phosphatase conjugated streptavidin.

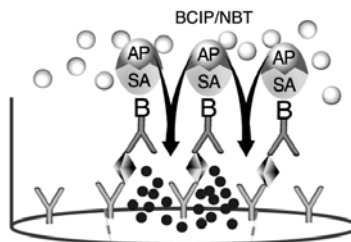

Add substrate and monitor the formation of colored spots.

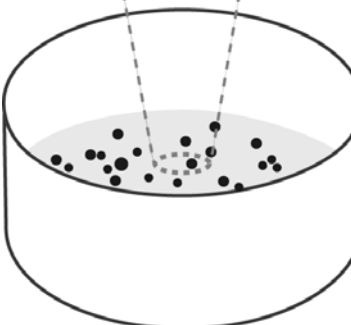

LIMITATIONS OF THE PROCEDURE

- FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES.
- The kit should not be used beyond the expiration date on the kit label.
- Any variation in pipetting and washing techniques, incubation time or temperature, or kit age can cause variation in density of spots, intensity of specific staining and background levels.

PRECAUTIONS

Some components of this kit contain sodium azide, which may react with lead and copper plumbing to form explosive metallic azides. Flush with large volumes of water during disposal.

Do not use reagents from this kit with components from other R&D Systems' ELISpot or ELISA kits and/or components manufactured by other vendors.

Do not remove the flexible plastic underdrain on the bottom of the microplate before or during incubation and development since it may damage the PVDF membrane filters. The underdrain cover may be removed only after completing the incubation with BCIP/NBT chromogen.

Although the toxicity of the chromogenic substrate BCIP/NBT is not currently known, wear gloves to avoid contact with skin. Follow local, state and federal regulations to dispose of used BCIP/NBT.

MATERIALS PROVIDED

Mouse TNF- α Microplate (Part 890900) - One 96-well PVDF-backed microplate coated with polyclonal antibody specific for mouse TNF- α .

Detection Antibody Concentrate (Part 890901) - 150 μ L of a 120X concentrated solution of biotinylated polyclonal antibody specific for mouse TNF- α with preservatives.

Streptavidin-AP Concentrate A (Part 895358) - 150 μ L of a 120X concentrated solution of Streptavidin conjugated to Alkaline Phosphatase with preservatives.

Dilution Buffer 1 (Part 895307) - 12 mL of a buffer for diluting Detection Antibody Concentrate with preservatives.

Dilution Buffer 2 (Part 895354) - 12 mL of a buffer for diluting Streptavidin-AP Concentrate A with preservatives.

10X Wash Buffer Concentrate (Part 895308) - 50 mL of a 10X concentrated solution of a buffered surfactant with preservative.

BCIP/NBT Chromogen (Part 895867) - 12 mL of a stabilized mixture of 5-Bromo-4-Chloro-3' Indolylphosphate p-Toluidine Salt (BCIP) and Nitro Blue Tetrazolium Chloride (NBT).

Mouse TNF- α Positive Control (Part 890902) - 2 ng of recombinant mouse TNF- α ; lyophilized.

STORAGE

Store the unopened kit at 2 - 8° C. Do not use beyond kit expiration date. This kit is validated for single use only. Results obtained with opened/reconstituted reagents at a later date may not be reliable.

OTHER SUPPLIES REQUIRED

- Pipettes and pipette tips
- Deionized or distilled water
- Multi-channel pipette, squirt bottle, manifold dispenser, or automated microplate washer
- 500 mL graduated cylinder
- 37° C CO₂ incubator
- Sterile culture media
- Dissection microscope or an automated ELISpot reader

TECHNICAL HINTS

- To minimize edge effect, place the microplate (bottom down) onto a piece of aluminum foil (about 4 x 6 inches). Add cells, cover the microplate with the lid and shape the foil around the edges of the microplate. The foil may be left on the microplate for the rest of the experimental procedure and removed after the BCIP/NBT has been washed off.
- Do not touch PVDF membrane filters with pipette tips when pipetting cells and reagents to avoid damage to the membrane.
- After completion of the experiment, do not dry the microplate at a temperature higher than 37° C since it may cause cracking of the PVDF membrane filters.
- The 96-well microplate provided in the kit is not sterile. However, due to the short incubation period and presence of antibiotics in the culture media, microbial contamination has not been a problem during the ELISpot procedure.
- The kit is designed for single use only. The layout of the assay should be carefully planned to maximize the use of the plate and reagents provided.
- The controls listed are recommended for each ELISpot experiment.
 - Positive Control - Use recombinant mouse TNF- α .
 - Unstimulated/Negative Control - Use the same number of unstimulated cells as stimulated cells.
 - Background Control - Use sterile culture media.
 - Detection Antibody Control - Substitute phosphate buffered saline for Detection Antibody.

REAGENT PREPARATION

Wash Buffer - If crystals have formed in the concentrate, warm to room temperature and mix gently until the crystals have completely dissolved. To prepare Wash Buffer, add 50 mL of 10X Wash Buffer Concentrate to 450 mL of deionized water and mix well.

Mouse TNF- α Positive Control - Reconstitute lyophilized mouse TNF- α with 250 μ L of culture medium that is used to incubate cells.

Detection Antibody - Tap or vortex the vial to release reagent collected in the cap. Transfer 100 μ L of Detection Antibody Concentrate into the vial labeled Dilution Buffer 1 and mix well. **For optimal performance, prepare Detection Antibody immediately before use.**

Streptavidin-AP - Tap or vortex the vial to release reagent collected in the cap. Transfer 100 μ L of Streptavidin-AP Concentrate A into the vial labeled Dilution Buffer 2 and mix well. **For optimal performance, prepare Streptavidin-AP immediately before use.**

SAMPLE PREPARATION

The types of effector and responder cells used, method of cell separation, mode of stimulation, and length of incubation are to be determined by each investigator. R&D Systems' cell selection products may be suitable for the purification of effector and responder cells. For a complete product listing of human, mouse, and rat cell selection products, see the R&D Systems' catalog or visit the web site at www.RnDSystems.com.

ASSAY PROCEDURE

Bring all reagents as needed to room temperature, except the Detection Antibody Concentrate and Dilution Buffer 1, which should remain at 2 - 8° C. All samples and controls should be assayed at least in duplicate. An Assay Record Template is provided at the back of this insert to record controls and samples assayed.

1. Fill all wells in the microplate with 200 μ L of sterile culture media and incubate for approximately 20 minutes at room temperature.
2. When cells are ready to be plated, aspirate the culture media from the wells. Immediately add 100 μ L of the appropriate cells or controls to each well (see Technical Hints for appropriate controls).
3. Incubate cells in a humidified 37° C CO₂ incubator. Optimal incubation time for each stimuli should be determined by the investigator. **Do not disturb the cells during the incubation period.**
4. Aspirate each well and wash, repeating the process three times for a total of four washes. Wash by filling each well with Wash Buffer (250 - 300 μ L) using a squirt bottle, multi-channel pipette, manifold dispenser, or autowasher. Complete removal of liquid at each step is essential to good performance. After the last wash, remove any remaining Wash Buffer by aspirating or decanting. Invert the plate and blot it against clean paper towels.
Note: Adjust the height of the prongs of the manifold dispenser or autowasher to prevent damage to the membranes.
5. Add 100 μ L of diluted Detection Antibody into each well and incubate at 2 - 8° C overnight.
6. Repeat step 4.
7. Add 100 μ L of diluted Streptavidin-AP into each well and incubate for 2 hours at room temperature.
8. Repeat step 4.
9. Add 100 μ L of BCIP/NBT Chromogen into each well and incubate for 1 hour at room temperature. **Protect from light.**
10. Discard the chromogen solution from the microplate and rinse the microplate with deionized water. Invert the microplate and tap to remove excess water. Remove the flexible plastic underdrain from the bottom of the microplate, wipe the bottom of the plate thoroughly with paper towels and dry completely either at room temperature (60 - 90 minutes) or 37° C (15 - 30 minutes).

CALCULATION OF RESULTS

The developed microplate can be analyzed by counting spots either manually using a dissection microscope or by using a specialized automated ELISpot reader. Specific spots are round and have a dark center with slightly fuzzy edges. Quantitation of results can be done, for example, by calculating the number of spot forming cells (SFC) per number of cells added into the well.

REPRODUCIBILITY DATA

Splenocytes from a C57BL mouse (1×10^6 cells/mL) were stimulated with 3 $\mu\text{g/mL}$ anti-mouse CD3 ϵ antibodies overnight at 37° C in a 5% CO₂ incubator. The sample was assayed in seven wells according to the procedure and analyzed with a dissection microscope.

Well	Number of Spots Counted
1	96
2	92
3	52
4	94
5	72
6	82
7	71

TROUBLESHOOTING GUIDE

Observation	Problem	Corrective Action
Following the incubation with BCIP/NBT chromogen and rinsing the microplate with deionized water, the dark-blue background color of filter membrane attenuates visualization and quantitation of spots.	Wet membrane	Microplates cannot be analyzed accurately until PVDF filter membranes are completely dry. Wait until membrane becomes dry, usually 15 - 30 minutes at 37° C or 60 - 90 minutes at room temperature.
The number of spots in the wells that contained the cells is high but their contrast as well as intensity of staining in the Positive Control wells is low.	Underdevelopment - perhaps as a result of using Streptavidin-AP and/or BCIP/NBT solutions that have not been adjusted to room temperature	Adjust the temperature of the reagents to room temperature before adding to the wells.
The number of spots in the wells that contained cells is lower than expected whereas Positive Control wells turned black-blue.	Cell stimulation problem	Ensure that reagents used to stimulate the cytokine release from the cells retained their biological activity. One way to check is to perform immunocytochemistry on fixed cells after stimulation.
	Too few cells added to the wells	Increase the number of cells added per well.
Following incubation with BCIP/NBT and drying the microplate, the density of the spots makes it difficult to quantify them.	Too many cells were added to the wells	Make dilutions of cells (<i>i.e.</i> , 1×10^6 , 5×10^5 , 1×10^5 , 5×10^4 , 1×10^4 cells per well) to determine the optimal number of cells that will result in formation of distinct spots.

REFERENCES

1. Kwon, B. *et al.* (1999) *Curr. Opin. Immunol.* **11**:340.
2. Idriss, H.T. and J.H. Naismith (2000) *Microsc. Res. Tech.* **50**:184.
3. Sedgwick, J.D. *et al.* (2000) *Immunol. Today* **21**:110.
4. Green, S. *et al.* (1976) *Proc. Natl. Acad. Sci. USA* **73**:381.
5. Mannel, D.N. *et al.* (1980) *Inf. Immun.* **30**:523.
6. Pennica, D. *et al.* (1984) *Nature* **312**:724.
7. Fransen, L. *et al.* (1985) *Nucl. Acids Res.* **13**:4417.
8. Kriegler, M. *et al.* (1988) *Cell* **53**:45.
9. Moss, M.L. *et al.* (1997) *Nature* **385**:733.
10. Marmenout, A. *et al.* (1985) *Eur. J. Biochem.* **152**:515.
11. Kwon, J. *et al.* (1993) *Gene* **132**:227.
12. Tchelingirian, J-L. *et al.* (1996) *J. Neurosci. Res.* **43**:99.
13. Lees, C.J. *et al.* (1999) *J. Interf. Cytokine Res.* **19**:1373.
14. Tashiro, M. *et al.* (1997) *J. Immunol.* **158**:2382.
15. Chensue, S.W. *et al.* (1988) *Am. J. Pathol.* **133**:564.
16. Medana, I.M. *et al.* (1997) *Am. J. Pathol.* **150**:1473.
17. Dahlen, E. *et al.* (1998) *J. Immunol.* **160**:3585.
18. Freyer, D. *et al.* (1999) *J. Immunol.* **163**:4308.
19. Gukovskaya, A.S. *et al.* (1997) *J. Clin. Invest.* **100**:1853.
20. Lewis, M. *et al.* (1991) *Proc. Natl. Acad. Sci. USA* **88**:2830.
21. Goodwin, R.G. *et al.* (1991) *Mol. Cell. Biol.* **11**:3020.
22. Baker, S.J. and E.P. Reddy (1998) *Oncogene* **17**:3261.
23. Vassalli, P. (1992) *Annu. Rev. Immunol.* **10**:411.
24. Paludan, S.R. (2000) *J. Leukoc. Biol.* **67**:18.
25. Czerkinsky, C.C. *et al.* (1983) *J. Immunol. Methods* **65**:109.
26. Sedgwick, J.D. and P.G. Holt (1983) *J. Immunol. Methods* **57**:301.
27. Czerkinsky, C.C. *et al.* (1984) *J. Immunol. Methods* **72**:489.
28. Helms, T. *et al.* (2000) *J. Immunol.* **164**:3723.

ASSAY RECORD TEMPLATE

This template may be used as a record of samples and controls run in an assay.

1								
2								
3								
4								
5								
6								
7								
8								
9								
10								
11								
12								
	A	B	C	D	E	F	G	H