

Monoclonal Antibody to CD44 - FITC

Alternate names:	CDw44, ECMR-III, Epican, Extracellular matrix receptor III, GP90 lymphocyte homing/adhesion receptor, HUTCH-I, Heparan sulfate proteoglycan, Hermes antigen, Hyaluronate receptor, LHR, MDU2, MDU3, MIC4, PGP-1, Phagocytic glycoprotein 1
Catalog No.:	CL110F
Quantity:	0.1 mg
Concentration:	0.1 mg/ml
Background:	<p>This antigen is expressed on most leukocytes (except a sub population of B cells) and increases upon activation. This antibody binds extracellularly to the standard (S) form on rat leukocytes, but it is not known if they bind to the N-terminal region. It has also been reported that the antibody may bind to melanoma cell lines that express CD44V (splice variant form).</p> <p>CD44 is expressed on most leukocytes except a sub population of B cells. Its expression is increased on T and B blasts.</p>
Uniprot ID:	P26051
NCBI:	10116
Host / Isotype:	Mouse / IgG2a
Clone:	OX-49
Immunogen:	T cell blasts.
Format:	State: Liquid purified IgG fraction Purification: Protein G affinity chromatography Buffer System: PBS with 0.02% sodium azide as preservative and EIA grade BSA as a stabilizer Label: FITC
Applications:	Flow cytometry (See protocol). Immunohistochemistry on Frozen Sections Immunohistochemistry on Paraffin Sections. Other applications not tested. Optimal dilutions are dependent on conditions and should be determined by the user.
Specificity:	This anti-Rat CD44 monoclonal antibody recognizes an epitope on both standard CD44 and its splice variant.
Species Reactivity:	Tested: Rat
Storage:	Store the antibody undiluted at 2-8°C for one month or (in aliquots) at -20°C for longer. Avoid repeated freezing and thawing. Shelf life: one year from despatch.

For research and in vitro use only. Not for diagnostic or therapeutic work.

Material Safety Datasheets are available at www.acris-antibodies.com or on request.

Antibody Hotline - Technical Questions - Antibody Location Service
Free Call: 0800-2274746 (Germany only) - www.acris-antibodies.com

- General References:**
1. Patterson, D.J., et al. 1987 Antigens of activated rat T lymphocytes including a molecule of 50,000 Mr detected only on CD4 positive T blasts. *Molec. Immunol.* 24(12): 1281-1290.
 2. Arch, R., et al. 1992. Participation in normal immune response of a metastases inducing splice variant of CD44. *Science.* 257:682-685.
 3. Wang, H., et al. 2001. Use of suppression subtractive hybridization for differential gene expression in stroke: discovery of CD44 gene expression and localization in permanent focal stroke in rats. *Stroke.* 32: 1020-1027.
 4. Jain, M., et al. 1996. Role of CD44 in the reaction of vascular smooth muscle cells to arterial wall injury. *J. Clin. Invest.* 97(3): 596-603.
 5. Lewington, A.J.P., et al. 2000. Expression of CD44 in kidney after acute ischemic injury in rats. *Am. J. Physiol.* 278: R247-R254.
 6. Foster, L.C., et al. 1998. Regulation of CD44 gene expression by the proinflammatory cytokine interleukin-1b in vascular smooth muscle cells. *J. Biol. Chem.* 273(32): 20341-20346.
 7. 4. *Stroke.* 32: 1020-1027.

Protocols:

FLOW CYTOMETRY ANALYSIS:

Method:

1. Prepare a cell suspension in media A. For cell preparations, deplete the red blood cell population with Rat cell separation medium.
2. Wash 2 times.
3. Resuspend the cells to a concentration of 2×10^7 cells/ml in media A. Add 50 μ l of this suspension to each tube (each tube will then contain 1×10^6 cells, representing 1 test).
4. To each tube, add 1.0-0.5 μ g* of CL110F or CL110FX.
5. Vortex the tubes to ensure thorough mixing of antibody and cells.
6. Incubate the tubes for 30 minutes at 4°C.
(It is recommended that the tubes are protected from light since most fluorochromes are light sensitive).
7. Wash 2 times at 4°C.
8. Resuspend the cell pellet in 50 μ l ice cold media B.
9. Transfer to suitable tubes for flow cytometric analysis containing 15 μ l of propidium iodide at 0.5 mg/ml in PBS. This stains dead cells by intercalating in DNA.

Media:

- A. Phosphate buffered saline (pH 7.2) + 5% normal serum of host species + sodium azide (100 μ l of 2M sodium azide in 100 mls).
- B. Phosphate buffered saline (pH 7.2) + 0.5% Bovine serum albumin + sodium azide (100 μ l of 2M sodium azide in 100 mls).

Results-Tissue Distribution:

Rat Strain: Wistar
Cell Concentration: 1×10^6 cells per tests
Antibody Concentration Used: 1.0 μ g/ 10^6 cells
Isotypic Control: FITC Mouse IgG2a.

Cell Source Percentage of cells stained above control:

Thymus: 99.6%
Spleen: 75.9%
Lymph Node: 96.8%

N.B. Appropriate control samples should always be included in any labelling studies.

* For optimal results in various applications, it is recommended that each investigator


For research and in vitro use only. Not for diagnostic or therapeutic work.

Material Safety Datasheets are available at www.acris-antibodies.com or on request.

Antibody Hotline - Technical Questions - Antibody Location Service
Free Call: 0800-2274746 (Germany only) - www.acris-antibodies.com

determine dilutions appropriate for individual use.

Pictures:

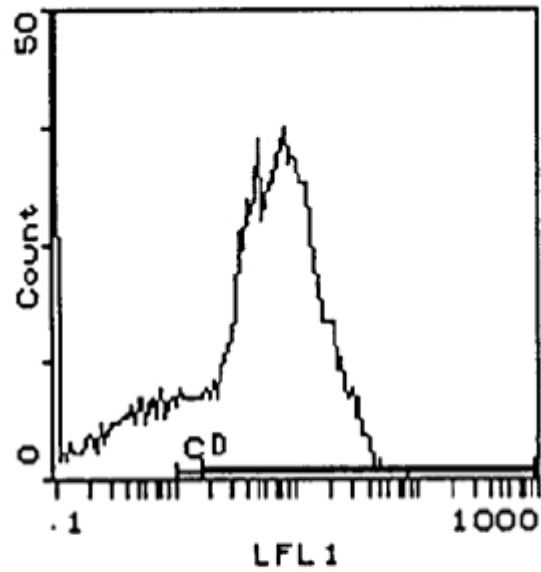

Cell Source: Spleen

Percentage of cells stained above control: 75.9%

For research and in vitro use only. Not for diagnostic or therapeutic work.

Material Safety Datasheets are available at www.acris-antibodies.com or on request.

Antibody Hotline - Technical Questions - Antibody Location Service
Free Call: 0800-2274746 (Germany only) - www.acris-antibodies.com