

Polyclonal Antibody to Collagen type II - Aff - Purified -

Alternate names:	Alpha-1 type II collagen, COL2A1
Catalog No.:	R1039X
Quantity:	0.5 mg
Concentration:	1.0 mg/ml (by UV absorbance at 280 nm)
Background:	Collagens are highly conserved throughout evolution and are characterized by an uninterrupted "Glycine-X-Y" triplet repeat that is a necessary part of the triple helical structure. For these reasons it is often extremely difficult to generate antibodies with specificities to collagens. The development of type specific antibodies is dependent on NON-DENATURED three-dimensional epitopes. For immunization collagens from human and bovine placenta and cartilage extensively purified by limited pepsin digestion and selective salt precipitation were used. This preparation results in a native conformation of the protein. Antibodies are isolated from rabbit antiserum and are extensively cross-adsorbed by immunoaffinity purification to produce 'type' specific antibodies. Greatly diminished reactivity and selectivity of these antibodies will result if denaturing and reducing conditions of SDS-PAGE and immunoblotting are used.
Uniprot ID:	P02458
NCBI:	NP_001835.3
GeneID:	1280
Host:	Rabbit
Immunogen:	Collagen type II purified from Human knee and Bovine nasal cartilage.
Format:	State: Liquid (sterile filtered) purified Ig fraction. Purification: Immunoaffinity Chromatography. Buffer System: 0.1 M Sodium Borate, 0.075M Sodium Chloride, 0.005 M EDTA, pH 8.0 Preservatives: 0.01% Sodium Azide Stabilizers: None
Applications:	Anti-Collagen antibodies have been used for Indirect trapping ELISA for quantitation of antigen in serum using a standard curve, for immunoprecipitation and for native (non-denaturing, non-dissociating) PAGE and Western blotting for highly sensitive qualitative analysis. <u>Recommended Dilutions:</u> ELISA: 1/5,000-1/50,000. Western Blot: 1/1,000-1/10,000. Immunoprecipitation: 1/100. Immunohistochemistry: 1/50-1/400. Note: This product was assayed by Immunoblot and found to be reactive against Collagen II at a dilution of 1/5,000 to 1/10,000.

For research and in vitro use only. Not for diagnostic or therapeutic work.

Material Safety Datasheets are available at www.acris-antibodies.com or on request.

Antibody Hotline - Technical Questions - Antibody Location Service
Free Call: 0800-2274746 (Germany only) - www.acris-antibodies.com

This product was also assayed against 1.0 µg of Collagen II in a standard Sandwich ELISA using Peroxidase conjugated affinity purified anti-Rabbit IgG [H&L] (Goat) and (ABTS (2,2'-azino-bis-[3-ethylbenthiazoline-6-sulfonic acid]) as a substrate for 30 minutes at room temperature. A working dilution of 1/4,000-1/8,000 of the stock concentration is suggested for this product.

Other applications not tested. Optimal dilutions are dependent on conditions and should be determined by the user.

Specificity:

This product has been prepared by Immunoaffinity Chromatography using immobilised antigens followed by extensive cross-adsorption against other collagens, Human serum proteins and non-collagen extracellular matrix proteins to remove any unwanted specificities.

Crossreactivity < 1% against other types of collagens was detected by ELISA against purified standards. Some class specific anti-collagens may be specific for three-dimensional epitopes which may result in diminished reactivity with denatured collagen or formalin-fixed, paraffin embedded tissues.

This antibody reacts with most mammalian Type II collagens and has negligible cross-reactivity with Type I, III, IV, V or VI collagens. Non-specific cross reaction of anti-collagen antibodies with other human serum proteins or non-collagen extracellular matrix proteins is negligible.

Species: Human, Mouse, Rat, Bovine, Sheep.

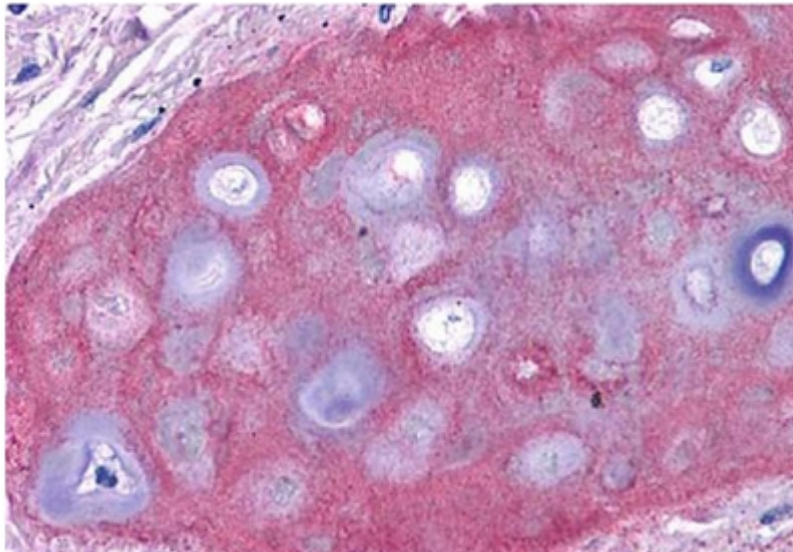
Other species not tested.

Storage:

Store the antibody undiluted at 2-8°C for one month or (in aliquots, mix with an equal volume of Glycerol) at -20°C for longer.

Avoid repeated freezing and thawing.

Shelf life: one year from despatch.

Pictures:

Collagen II antibody Cat.-No R10369 (Lot 26014, 1/400, 45 min RT) showed moderate staining of in FFPE Human bronchiolar cartilage (shown in image). Though not shown, faint to moderate staining of tonsillar squamous epithelium, prostatic stroma, breast, colon, placenta, and dermal connective tissues was also observed. All other tissues, including brain, breast epithelium, colon epithelium, heart, intestine, kidney, liver, lung, skeletal muscle, pancreas, spleen, testis, thymus, thyroid, and uterus were negative for staining. Slides were steamed in 0.01M sodium citrate buffer, pH 6.0 at 99-100°C - 20 minutes for antigen retrieval. Image provided courtesy of LifeSpan Biosciences, Seattle, WA

For research and in vitro use only. Not for diagnostic or therapeutic work.

Material Safety Datasheets are available at www.acris-antibodies.com or on request.

Antibody Hotline - Technical Questions - Antibody Location Service
Free Call: 0800-2274746 (Germany only) - www.acris-antibodies.com