

MONKEY MYOGLOBIN ELISA

Life Diagnostics, Inc., Catalog Number: MYO-3

ELISA for the Determination of Myoglobin in Monkey Serum, Plasma & Urine¹

STORAGE

Store standard at -20°C
STORE REMAINDER OF KIT AT 2 - 8°C

INTRODUCTION

Myoglobin is a heme protein found in both cardiac and skeletal muscle. Following myocardial necrosis associated with myocardial infarction (MI), myoglobin is one of the first markers to rise above normal levels. Studies with human subjects have shown that myoglobin increases measurably above baseline within 2-4 hours post-infarct, peaking at 9-12 hours, and returning to baseline within 24-36 hours. In the absence of skeletal muscle trauma or other factors associated with a noncardiac related increase in circulating myoglobin, myoglobin may be used as a marker for MI. Similarly, in the absence of cardiac damage, myoglobin may be used as a marker of skeletal muscle injury.

PRINCIPLE OF THE ASSAY

The myoglobin test kit is based on the principle of a solid phase enzyme-linked immunosorbent assay (ELISA). A monoclonal anti-myoglobin antibody is used for solid phase immobilization (on the microtiter wells) and a polyclonal anti-myoglobin antibody conjugated to horseradish peroxidase (HRP) is used for detection. The test sample is allowed to react simultaneously with the two antibodies, and the myoglobin molecules are thus sandwiched between the solid phase and enzyme-linked antibodies. After one hour incubation at room temperature, the wells are washed to remove unbound labeled antibodies. TMB (tetramethylbenzidine) Reagent is added and incubated for 20 minutes, resulting in the development of a blue color. The color development is stopped with the addition of Stop Solution, changing the color to yellow. The concentration of myoglobin in the test sample is proportional to the absorbance measured spectrophotometrically at 450 nm.

MATERIALS AND COMPONENTS

Materials provided with the kit:

- Anti-myoglobin coated microtiter wells, 96 wells
- Monkey Myoglobin Standard (100 µl, 50 µg/ml)²
- Diluent, 12 ml
- Enzyme Conjugate Reagent, 11 ml
- 20x Wash Solution, 50 ml
- TMB Reagent (One-Step), 11 ml

¹ A matrix effect may be observed with urine samples, resulting in slight variances in absorbance values relative to dilution factor. It is therefore recommended that, wherever possible, all urine samples within a particular study be similarly diluted (with the provided diluent) prior to testing, thereby ensuring accurate determination of relative myoglobin levels.

² Due to international import/export restrictions of monkey derived products, the myoglobin standard provided with this kit is of non monkey origin. The standard curve obtained with this material is identical to that obtained with monkey myoglobin.

- Stop Solution (1N HCl), 11 ml
- Materials required but not provided:**

- Precision pipettes
- Disposable pipette tips
- Distilled or deionized water
- Vortex mixer
- Absorbent paper or paper towels
- Graph paper (PC graphing software is optional but recommended)
- Plate shaker
- Microtiter plate reader with an optical density range of 0-4 OD at 450 nm wavelength

WASH SOLUTION PREPARATION

The wash solution is provided as a 20x stock. Prior to use, dilute the contents of the bottle (50 ml) with 950 ml of distilled or deionized water.

STANDARD PREPARATION

1. Label 8 polypropylene tubes as 125, 62.5, 31.25, 15.63, 7.81, 3.91, 1.95 and 0 ng/ml.
2. Pipette 997.5 µl of diluent into the tube labeled 125 ng/ml.
3. Pipette 100 µl of diluent into the tubes labeled 62.5, 31.25, 15.63, 7.81, 3.91, 1.95 and 0 ng/ml.
4. Briefly centrifuge or flick the myoglobin stock tube to ensure that the liquid contents are at the bottom of the tube.
5. Dilute 2.5 µl of the 50 µg/ml myoglobin stock into the 997.5 µl of diluent in the tube labeled 125 ng/ml. This provides a 125 ng/ml solution of myoglobin.
6. Prepare a 62.5 ng/ml stock by diluting and mixing 100 µl of the 125 ng/ml stock with 100 µl of diluent in the tube labeled 62.5 ng/ml. Similarly prepare the 31.25, 15.63, 7.81, 3.91 and 1.95 ng/ml stocks by serial dilution.

SAMPLE COLLECTION

Serum, plasma and urine should be collected using standard techniques. Remove serum or plasma from the coagulated or packed cells within 60 minutes after collection. Plasma samples may be collected into tubes containing EDTA. Samples that cannot be assayed within 3 hours of collection should be frozen at -20°C or lower. Samples should not be repeatedly frozen and thawed prior to testing.

SAMPLE PREPARATION

Samples may be tested undiluted or after dilution with diluent. The dilution factor should be determined empirically. *Only 20 µl of sample is required per assay (2 x 20 µl, if samples are to be tested in duplicate).*

ASSAY PROCEDURE

1. Ensure that all reagents are at room temperature.
2. Secure the desired number of coated wells in the holder.
3. Dispense 100 µl of Enzyme Conjugate Reagent into each well.
4. Dispense 20 µl of myoglobin standards and samples (in duplicate) into the appropriate wells.

5. Incubate at room temperature (18-25°C) on a plate shaker for one hour. Mix Gently (~100-150 rpm)
6. Remove the incubation mixture either with a plate washer or by flicking plate contents into a waste container.
7. Wash and empty the microtiter wells 5 times with 1x wash solution. This may be performed using either a plate washer (400 µl/well) or a squirt bottle. The entire wash procedure should be performed as quickly as possible.
8. Strike the wells sharply onto absorbent paper or paper towels to remove all residual water drops.
9. Dispense 100 µl of TMB Reagent solution into each well. Gently mix for 5 seconds.
10. Incubate on a plate shaker at room temperature for 20 minutes. Mix Gently.
11. Stop the reaction by adding 100 µl of Stop Solution to each well.
12. Gently mix for 30 seconds. It is important to make sure that all the blue color changes to yellow color completely.
13. Read absorbance at 450 nm with a microtiter well reader within 5 minutes.


CALCULATION OF RESULTS

1. Calculate the mean absorbance value for each set of reference standards and samples.
2. Construct a standard curve by plotting the mean absorbance obtained for each reference standard against its concentration in ng/ml on graph paper, with absorbance values on the vertical or Y axis and concentrations on the horizontal or X axis.
3. Use the mean absorbance values for each sample to determine the corresponding concentration of myoglobin in ng/ml from the standard curve.
4. Multiply the derived value by the appropriate dilution factor if the test samples were diluted.
5. Graphing software, if available, should be used.

TYPICAL STANDARD CURVE

Results of a typical standard run with optical density readings at 450 nm on the Y axis against myoglobin concentrations on the X axis are illustrated below. This standard curve is for illustrative purposes only and should not be used to calculate unknowns. A standard curve should be run for each assay.

Myoglobin (ng/ml)	Absorbance (450 nm)
125	2.320
62.5	1.358
31.25	0.738
15.63	0.441
7.81	0.276
3.91	0.190
1.98	0.158
0	0.111

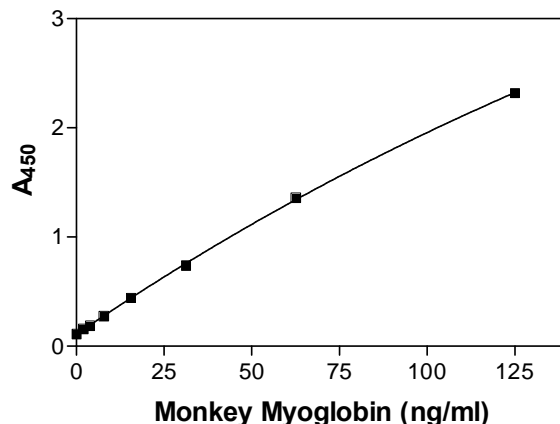

PROCEDURAL NOTES

1. Standards should be prepared immediately prior to use and should be used within 30 minutes of preparation.
2. Pipetting of conjugate, standards and samples into the microtiter plate should be completed within 10 minutes.
3. We recommend that standards and samples be run in duplicate.
4. It is recommended that the wells be read within 5 minutes following addition of Stop Solution.
5. Reliable and reproducible results will be obtained when the assay procedure is carried out with a complete understanding of the package insert instructions.
6. The wash procedure is critical. Insufficient washing will result in poor precision and falsely elevated absorbance readings

Rev 041015NC

For technical assistance please email us at techsupport@lifediagnosics.com