



Rhodopsin / RHO Mouse anti-Human Monoclonal (B630) Antibody - LS-B10408 - LSBio		
CatalogID:	LS-B10408	
Validation:	This antibody replaces catalog number LS-C204581. It has been validated for use in the following assays: IHC-P.	
Target:	rhodopsin	
Synonyms:	RHO Antibody, CSNBAD1 Antibody, OPN2 Antibody, Rhodopsin Antibody, RP4 Antibody, Opsin 2, rod pigment Antibody, Opsin-2 Antibody	
Family / Subfamily:	GPCR / Opsin	
Host	RHO antibody was produced in Mouse	
Clonality:	Monoclonal	
Isotype:	IgG1,k	
Clone Name:	B630	
Immunogen Species:	Rhodopsin / RHO antibody was raised against Human	
Reactivity:	Human, Mouse, Rat, Bovine, Pig	
Purification:	Protein G purified	
Presentation:	PBS, 10 mM sodium azide.	
Recommended Storage:	+4°C or -20°C, Avoid repeated freezing and thawing.	
Usage Summary:	Try at dilutions of ~1:1000 for immunofluorescence. For western blots try at 1:5000. A suitable control tissue is retinal extracts. Rhodopsin run at 35 kDa on SDS-PAGE gels.	
Uses:	IHC - Paraffin (10 μg/ml), ICC (1:1000), Immunofluorescence (1:1000), Western blot (1:5000) (Optimal dilution to be determined by the researcher)	
Size:	50 µl	

Immunohistochemistry Image:

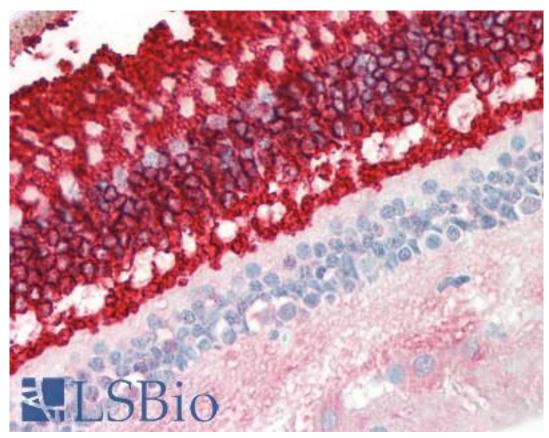


Human Retina: Formalin-Fixed, Paraffin-Embedded (FFPE)

Immunocytochemistry Image:

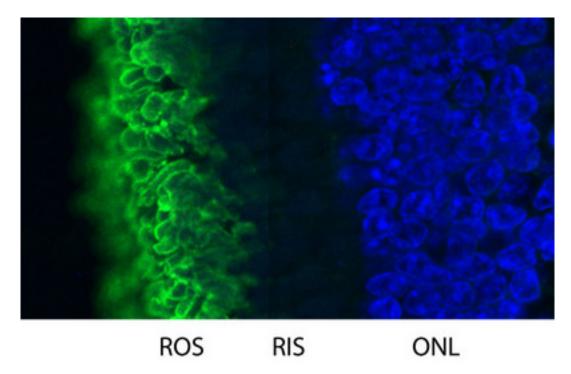

High magnification confocal image of pig retinal section stained with Rhodopsin / RHO antibody (Green). Rhodopsin is most abundant in the rod outer segments (ROS) of retina, clearly localized in rod membranes. The rod inner segments (RIS) and rod nuclei in the outer nuclear layer (ONL) are also seen in this image. Nuclear DNA was stained with DAPI (blue).

Western Blot Image:

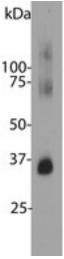

Blot of bovine retinal extracts probed with Rhodopsin / RHO antibody. The antibody stains a band corresponding to retinal rhodopsin at about 35kDa. Bands about 70 kDa and 140 kDa are aggregated forms of rhodopsin. Note, due to the highly hydrophobic nature of rhodopsin, it important to avoid boiling samples containing this protein it in SDS-PAGE sample buffer, as this will result in even more extensive aggregation of the rhodopsin protein and appearance of more of this high molecular weight material.

Requested From:	Japan
-----------------	-------

Laboratory Reagent For In Vitro Research Use Only
Not for resale without prior written consent from LifeSpan BioSciences, Inc.
Created on 9/23/2014
© 2014 LifeSpan BioSciences