


# **RayBio®Phospho- Akt (Ser473) ELISA Kit**

**For Measuring Phosphorylated Akt (Ser473)  
in Human, Mouse and Rat Cell Lysates**

**User Manual  
(Revised Mar 1, 2012)**

**RayBio® Akt (Ser473) ELISA Kit  
Protocol**

(Cat#: PEL-Akt-S473)

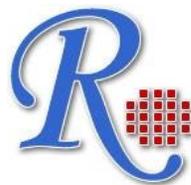


**RayBiotech, Inc.**

**We Provide You With Excellent  
Protein Array System And Service**

**Tel:(Toll Free)1-888-494-8555 or 770-729-2992; Fax:770-206-2393;**

**Web: [www.raybiotech.com](http://www.raybiotech.com) Email: [info@raybiotech.com](mailto:info@raybiotech.com)**

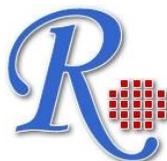

RayBiotech, Inc.

---

**RayBio® Phospho-Akt (Ser473)  
ELISA Kit Protocol**

**TABLE OF CONTENTS**

I.	Introduction.....	2
II.	Material Provided.....	2
III.	Storage.....	3
IV.	Additional Materials Required.....	4
V.	Sample Preparation.....	4
VI.	Reagent Preparation.....	5
VII.	Assay Procedure.....	7
VIII.	Assay Procedure Summary.....	8
IX.	Typical Data.....	9
	i. Positive Control.....	9
	ii. Recombinant Human PDGF Stimulation of NIH3T3 Cell Lines.....	10
	iii. Sensitivity.....	11
X.	References.....	11
XI.	Troubleshooting Guide.....	12

## **I. INTRODUCTION**

RayBio®Phospho-Akt (Ser473) ELISA (Enzyme-Linked Immunosorbent Assay) kit is a very rapid, convenient and sensitive assay kit that can monitor the activation or function of important biological pathways in human and mouse cell lysates. By determining phosphorylated Akt protein in your experimental model system, you can verify pathway activation in your cell lysates. You can simultaneously measure numerous different cell lysates without spending excess time and effort in performing a Western Blotting analysis.

This Sandwich ELISA kit is an in vitro enzyme-linked immunosorbent assay for the measurement of human, mouse and rat phospho-Akt (Ser473). An anti-pan Akt antibody has been coated onto a 96-well plate. Samples are pipetted into the wells and Akt present in a sample is bound to the wells by the immobilized antibody. The wells are washed and anti-Akt (Ser473) antibody is used to detect phosphorylated Akt (Ser473). After washing away unbound antibody, HRP-conjugated anti-rabbit IgG is pipetted to the wells. The wells are again washed, a TMB substrate solution is added to the wells and color develops in proportion to the amount of Akt (Ser473) bound. The Stop Solution changes the color from blue to yellow, and the intensity of the color is measured at 450 nm.

## **II. MATERIAL PROVIDED**

1. Akt Microplate (Item A): 96 wells (12 strips x 8 wells) coated with anti-pan Akt antibody.

2. Wash Buffer Concentrate (20x) (Item B): 25 ml of 20x concentrated solution
3. Assay Diluent (Item E2): 15 ml of 5x concentrated buffer. For diluting cell lysate sample, detection antibody (Item C) and HRP-conjugated anti-rabbit IgG Concentrate (Item D-1).
4. Detection Antibody Akt (Ser473) (Item C): 2 vial of rabbit anti-Akt (Ser473) (each vial is enough to assay half microplate).
5. HRP-conjugated Anti-rabbit IgG (Item D-1), 25  $\mu$ l of 500x HRP-conjugated anti-rabbit IgG concentrate.
6. TMB One-Step Substrate Reagent (Item H): 12 ml of 3,3',5,5'-tetramethylbenzidine (TMB) in buffered solution.
7. Stop Solution (Item I): 8 ml of 0.2 M sulfuric acid.
8. Cell Lysate Buffer (Item J): 5 ml 2x cell lysis buffer (not including protease and phosphatase inhibitors).
9. Positive Control A431S003-1 (Item K): 1 vial of lyophilized powder from A431 cell lysate.

### **III. STORAGE**

Upon receipt, the kit should be stored at  $-20^{\circ}\text{C}$ . Please use within 6 months from the date of shipment. After initial use, Wash Buffer Concentrate (Item B), Assay Diluent (Item E2), TMB One-Step Substrate Reagent (Item H), Stop Solution (Item I) and Cell Lysate Buffer (Item J) should be stored at  $4^{\circ}\text{C}$  to avoid repeated freeze-thaw cycles. Return unused wells to the pouch containing desiccant pack, reseal along entire edge and store at  $-20^{\circ}\text{C}$ . Item D-1 store at  $2-8^{\circ}\text{C}$  for up to one month (store at  $-20^{\circ}\text{C}$  for up to 6 months, avoid repeated freeze-thaw cycles). Reconstituted Positive Control (Item K) should be stored at  $-70^{\circ}\text{C}$ .

## **IV. ADDITIONAL MATERIALS REQUIRED**

- 1 Microplate reader capable of measuring absorbance at 450 nm.
- 2 Protease and Phosphatase inhibitors.
- 3 Shaker.
- 4 Precision pipettes to deliver 2  $\mu$ l to 1 ml volumes.
- 5 Adjustable 1-25 ml pipettes for reagent preparation.
- 6 100 ml and 1 liter graduated cylinders.
- 7 Distilled or deionized water.
- 8 Tubes to prepare sample dilutions.

## **V. SAMPLE PREPARATION**

Cell lysates - Rinse cells with PBS, making sure to remove any remaining PBS before adding the Cell Lysate Buffer. Solubilize cells at  $4 \times 10^7$  cells/ml in 1x Cell Lysate Buffer (we recommend adding protease and phosphatase inhibitors to Cell Lysate Buffer prior to sample preparation). Pipette up and down to resuspend and incubate the lysates with shaking at 2 - 8° C for 30 minutes.

Microcentrifuge at 13,000 rpm for 10 minutes at 2 - 8° C, and transfer the supernates into a clean test tube. Lysates should be used immediately or aliquoted and stored at -70 °C. Avoid repeated freeze-thaw cycles. Thawed lysates should be kept on ice prior to use.

For the initial experiment, we recommend to do a serial dilution testing such as 5-fold and 50-fold dilution for your cell lysates with 1x Assay Diluent (Item E2) before use.

*Note: The fold dilution of sample used depends on the abundance of phosphorylated proteins and should be determined empirically. More of the sample can be used if signals are too weak. If signals are too strong, the sample can be diluted further.*


Cell Lysate Buffer should be diluted 2-fold with deionized or distilled water before use (recommend to add protease and phosphatase inhibitors).

## **VI. REAGENT PREPARATION**

1. Bring all reagents and samples to room temperature (18 - 25°C) before use.
2. Item E2, Assay Diluent should be diluted 5-fold with deionized or distilled water before use.
3. Preparation of Positive Control: Briefly spin the Positive Control vial of Item K. Add 500 µl 1x Assay Diluent (Item E2, Assay Diluent should be diluted 5-fold with deionized or distilled water before use) into Item K vial to prepare Positive Control (1) Solution (See i. Positive Control of part IX. TYPICAL DATA for a typical result in page 9). **Dissolve the powder thoroughly by a gentle mix** (it can be removed by centrifuge if any precipitate in the solution is found). Pipette 300 µl 1x Assay Diluent into each tube. Use the Positive Control (1) to produce a dilution series (shown below). Mix each tube thoroughly before the next transfer. 1x Assay Diluent serves as the background.

Positive Control powder +  
500  $\mu$ l 1x Assay Diluent

150 $\mu$  150 150  $\mu$ l

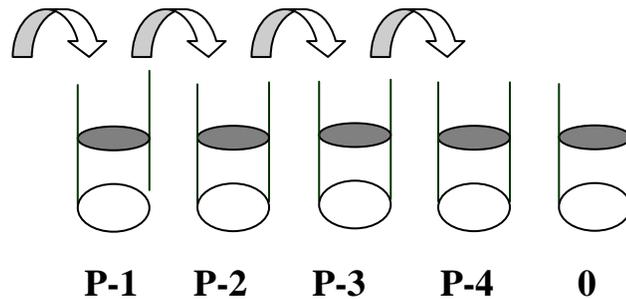

4. If the Wash Concentrate (20x) (Item B) contains visible crystals, warm to room temperature and mix gently until dissolved. Dilute 20 ml of Wash Buffer Concentrate into deionized or distilled water to yield 400 ml of 1x Wash Buffer.
5. Briefly spin the detection antibody (Item C) before use. Add 100  $\mu$ l of 1x Assay Diluent into the vial to prepare a detection antibody concentrate. Pipette up and down to mix gently (the concentrate can be stored at 4°C for 5 days or at – 80°C for one month). The rabbit anti-Akt (Ser473) antibody should be diluted 55-fold with 1x Assay Diuent and used in step 4 of Part VII Assay Procedure.
6. Briefly spin the HRP-conjugated anti-rabbit IgG (Item D) before use. Pipette up and down to mix gently. HRP-conjugated anti-rabbit IgG concentrate should be diluted 500-fold with 1x Assay Diuent.

*For example: Briefly spin the vial (ItemD) and pipette up and down to mix gently. Add 10  $\mu$ l of HRP-conjugated anti-*

*rabbit IgG concentrate into a tube with 5 ml 1x Assay Diluent to prepare a 500-fold diluted HRP-conjugated anti-rabbit IgG solution.*

7. Cell Lysate Buffer should be diluted 2-folds with deionized or distilled water before use (recommend to add protease and phosphatase inhibitors).

## **VII. ASSAY PROCEDURE:**

1. Bring all reagents to room temperature (18 - 25°C) before use. It is recommended that all samples or Positive Control should be run at least in duplicate.
2. Add 100 µl of each sample or positive control into appropriate wells. Cover well with plate holder and incubate for 2.5 hours at room temperature or over night at 4°C with shaking.
3. Discard the solution and wash 4 times with 1x Wash Solution. Wash by filling each well with Wash Buffer (300 µl) using a multi-channel pipette or autowasher. Complete removal of liquid at each step is essential to good performance. After the last wash, remove any remaining Wash Buffer by aspirating or decanting. Invert the plate and blot it against clean paper towels.
4. Add 100 µl of prepared 1x detection antibody anti-Akt (Ser473) (Reagent Preparation step 5) to each well. Incubate for 1 hour at room temperature with shaking.

5. Discard the solution. Repeat the wash as in step 3.
6. Add 100  $\mu$ l of prepared 1x HRP-conjugated anti-rabbit IgG (see Reagent Preparation step 6) to each well. Incubate for 1 hour at room temperature with shaking.
7. Discard the solution. Repeat the wash as in step 3.
8. Add 100  $\mu$ l of TMB One-Step Substrate Reagent (Item H) to each well. Incubate for 30 minutes at room temperature in the dark with shaking.
9. Add 50  $\mu$ l of Stop Solution (Item I) to each well. Read at 450 nm immediately.

## **VIII. ASSAY PROCEDURE SUMMARY**

1. Prepare all reagents, samples and standards as instructed.  

2. Add 100  $\mu$ l sample or positive control to each well.  
Incubate 2.5 hours at room temperature or over night at 4°C.  

3. Add 100  $\mu$ l prepared primary antibody to each well.  
Incubate 1.0 hours at room temperature.  

4. Add 100  $\mu$ l prepared 1X HRP-Streptavidin solution.  
Incubate 1 hour at room temperature.  


5. Add 100  $\mu$ l TMB One-Step Substrate Reagent to each well.  
Incubate 30 minutes at room temperature.




6. Add 50  $\mu$ l Stop Solution to each well.  
Read at 450 nm immediately.

## IX. TYPICAL DATA

ELISA data analysis: Average the duplicate readings for each sample or positive.

### i. Positive Control


A431 cells were treated with recombinant human EGF at 37°C for 20 min. Solubilize cells at  $4 \times 10^7$  cells/ml in Cell Lysate Buffer. Serial dilutions of lysates were analyzed in this ELISA. Please see step 3 of Part VI Reagent Preparation for detail.

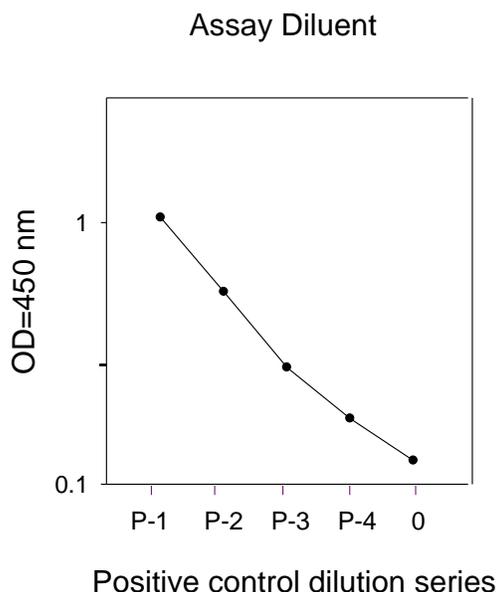


## ii. Recombinant Human PDGF Stimulation of NIH3T3 Cell Lines

NIH3T3 cells were treated or untreated with recombinant human PDGF for 10 min. Cell lysates were analyzed using this phosphoELISA and Western Blot.

### A). ELISA

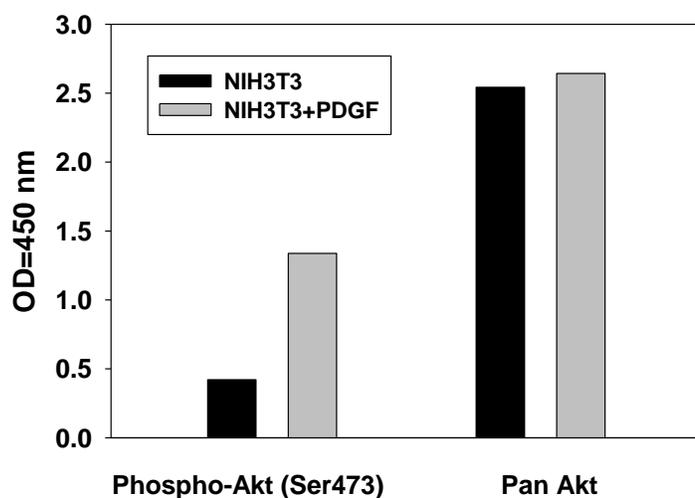


### B). Western-Blot Analysis

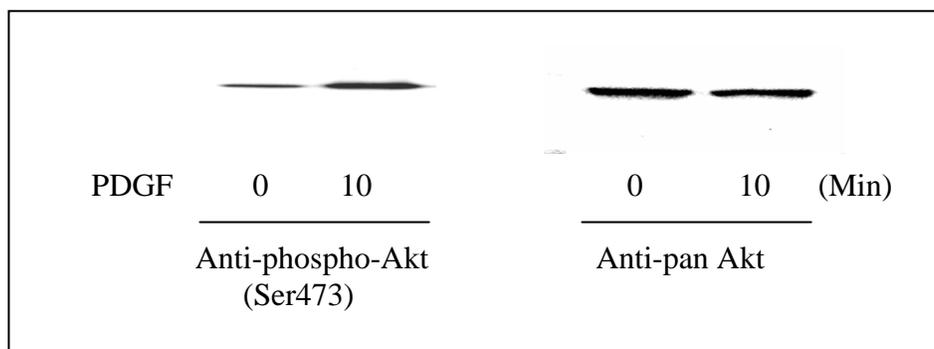

### iii. SENSITIVITY

The NIH3T3 cells were treated with recombinant human PDGF for 10 minutes to induce phosphorylation of Akt. Serial dilutions of lysates were analyzed in this ELISA and by Western blot. Immunoblots were incubated with anti-phospho-Akt (Ser473).

#### A) ELISA

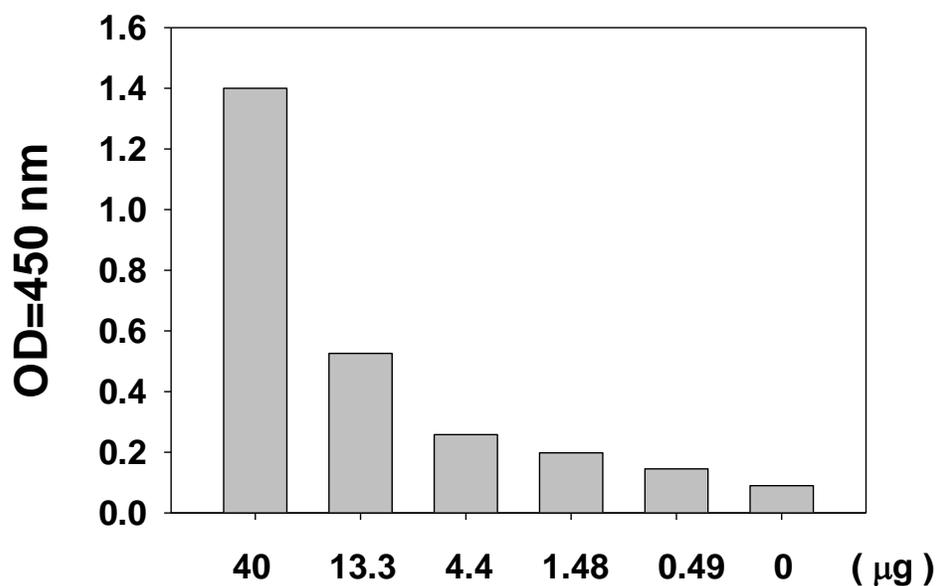

### X. REFERENCES:

1. Hajduch, E. et al. (2001) *FEBS Lett*, 492, 199-203.
2. Burgering, B.M. and Coffey, P.J. (1995) *Nature*, 376, 599-602.
3. Franke, T.F. et al. (1997) *Cell*, 88, 435-7.

## XI. TROUBLESHOOTING GUIDE

<b>Problem</b>	<b>Cause</b>	<b>Solution</b>
1. Sample signals: a. Too low  b. Too high	a. Sample concentration is too low  b. Sample concentration is too high	a. Increasing sample concentration  b. Reducing sample concentration
2. Large CV	a. Inaccurate pipetting	a. Check pipettes
3. High background	a. Plate is insufficiently washed  b. Contaminated wash buffer	a. Review the manual for proper washing. If using an automated plate washer, check that all ports are unobstructed.  b. Make fresh wash buffer
4. Positive Control: Low signal	a. Improper storage of the ELISA kit  b. Stop solution  c. Improper primary or secondary antibody dilution	a. Upon receipt, the kit should be stored at $-20^{\circ}\text{C}$ . Store the positive control at $-70^{\circ}\text{C}$ after reconstitution.  b. Stop solution should be added to each well before measurement and read OD immediately.  c. Ensure correct dilution

RayBio® ELISA kits:

Over 200 ELISA kits, custom ELISA kit choose from over 500 list visit [www.raybiotech.com](http://www.raybiotech.com) for details.

RayBiotech, Inc., the protein array pioneer company, strives to research and develop new products to meet demands of the biomedical community. RayBio's patent-pending technology allows detection of over 180 cytokines, chemokines and other proteins in a single experiment. Our format is simple, sensitive, reliable and cost effective. Products include: Cytokine Arrays, Chemokine Arrays, ELISA kits, Phosphotyrosine kits, Recombinant Proteins, Antibodies, and custom services.

### **Antibody Array**

**Cytokine Antibody Array: Simultaneous detection up to 200 proteins (cytokine, chemokine, growth factor, adipokine, angiogenic factor, protease) in one experiment**

### **Phosphorylation Antibody Array**

- RTK antibody array
- STAT1 phosphorylation antibody arrays

**Label based antibody array: Simultaneous detection more than 500 proteins in one experiment**

**Quantibody Array: Quantitative measurement of multiple protein levels**

### **Protein Array**

### **ELISA**

### **Cell-Based Phosphorylation ELISA**

### **Tissue MicroArray**

**Protein: Cytokine, Chemokine, Adipokine, Angiogenic factor, Virus, bacteria and infectious disease protein, hormone, Enzyme, other**

### **Peptide**

**Antibody: Cytokine, Adipokine, Angiogenic factor, Signal transduction, Transcription factor, Receptor, Adhesion molecule, Virus, bacteria and other infectious agents, Secondary antibody, Tag antibody, Immunoglobulin, Hormone, Cell surface, Protease, other**

**Antibody array, Protein array, Peptide array, ELISA, Phosphorylation assay  
Tissue array**

**Assay service:** just simply send your samples and get data in 1 to 2 weeks.

Antibody array, Protein array, ELISA, Quantibody array

**Antibody production:** highest quality with very competitive price

Monoclonal antibody, Recombinant antibody, Polyclonal antibody, Phase display, Antibody engineering, Antibody conjugation

**Recombinant protein production**

**Assay development**


**Array printing**

Contact and non-contact arrayers. All kinds of substrates of your choice including glass slides, membranes and plates.

**Note:**

**Note:**

This product is for research use only.

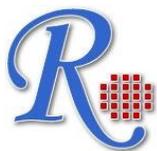

©2004 *RayBiotech, Inc.*