DuoSet® IC

Human Total XIAP

Catalog Number DYC822-2

DYC822-5

DYC822E

For the development of sandwich ELISAs to measure XIAP in cell lysates.

This package insert must be read in its entirety before using this product.

FOR RESEARCH USE ONLY.
NOT FOR USE IN DIAGNOSTIC PROCEDURES.

TABLE OF CONTENTS

Contents	Page
PRINCIPLE OF THE ASSAY	2
MATERIALS PROVIDED	2
OTHER MATERIALS REQUIRED	3
SOLUTIONS REQUIRED	3
REAGENT PREPARATION	4
PREPARATION OF SAMPLES	4
PRECAUTION	5
TECHNICAL HINTS AND LIMITATIONS	5
GENERAL ELISA PROTOCOL	6
CALCULATION OF RESULTS	7
TYPICAL DATA	7
CALIBRATION	8
SPECIFICITY	8
QUANTIFICATION	9
PLATE LAYOUT	10

MANUFACTURED AND DISTRIBUTED BY:

R&D Systems, Inc. TELEPHONE: (800) 343-7475
614 McKinley Place NE (612) 379-2956
Minneapolis, MN 55413 FAX: (612) 656-4400

United States of America E-MAIL: info@RnDSystems.com

DISTRIBUTED BY:

R&D Systems Europe, Ltd.

 19 Barton Lane
 TELEPHONE:
 +44 (0)1235 529449

 Abingdon Science Park
 FAX:
 +44 (0)1235 533420

 Abingdon, OX14 3NB
 E-MAIL:
 info@RnDSystems.co.uk

United Kingdom

R&D Systems China Co. Ltd.

24A1 Hua Min Empire Plaza TELEPHONE: +86 (21) 52380373 726 West Yan An Road FAX: +86 (21) 52371001

Shanghai PRC 200050 E-MAIL: info@RnDSystemsChina.com.cn

PRINCIPLE OF THE ASSAY

This DuoSet[®] IC ELISA contains the basic components required for the development of sandwich ELISAs to measure XIAP in cell lysates. An immobilized capture antibody binds XIAP present in samples or standards. After washing away unbound material, a biotinylated detection antibody specific for XIAP is used to detect the protein, utilizing a standard Streptavidin-HRP format.

MATERIALS PROVIDED

Store the unopened kit at 2-8° C. Do not use past kit expiration date.

		Vials Provided		
Description	Part #	Storage Conditions	Cat. # DYC822-2	Cat. # DYC822-5
Total XIAP Capture Antibody	840758	2-8° C	1	2
Total XIAP Detection Antibody	840759	2-8° C	1	2
Total XIAP Standard	840760	2-8° C	3	5
Streptavidin-HRP	890803	2-8° C	1	1

DYC822-2 contains sufficient materials to run ELISAs on at least two 96 well plates.*

DYC822-5 contains sufficient materials to run ELISAs on at least five 96 well plates.*

This kit is also available in an Economy Pack (R&D Systems, Catalog # DYC822E). Economy Packs contain sufficient materials to run ELISAs on 15 microplates. Specific vial counts of each component may vary. Please refer to the literature accompanying your order for specific vial counts.

*Provided the following conditions are met:

- The reagents are prepared as described in this package insert.
- The assay is run as described in the General ELISA Protocol on page 6.
- The recommended microplates, buffers, diluents, substrates, and solutions are used.

OTHER MATERIALS REQUIRED

- Aprotinin (Sigma # A6279)
- Leupeptin (Sigma # L8511)
- Pepstatin (Sigma # P4265)
- Phenylmethylsulfonylfluoride (PMSF) (Sigma # P7626)
- Sodium Azide (NaN₃) (Sigma # S2002)
- Triton® X-100 (Sigma # T9284)
- Pipettes and pipette tips
- Deionized or distilled water
- 96 well microplates [Costar EIA Plates (Catalog # 2592 or R&D Systems Catalog # DY990) are suggested]
- Plate sealers (R&D Systems, Catalog # DY992)
- Squirt bottle, manifold dispenser, or automated microplate washer.

SOLUTIONS REQUIRED

PBS - 137 mM NaCl, 2.7 mM KCl, 8.1 mM Na₂HPO₄, 1.5 mM KH₂PO₄, pH 7.2-7.4, 0.2 μ m filtered.

Wash Buffer - 0.05% Tween $^{\mathbb{R}}$ 20 in PBS, pH 7.2-7.4 (R&D Systems, Catalog # WA126).

Block Buffer - 1% BSA*, 0.05% NaN₃, in PBS, pH 7.2-7.4.

IC Diluent #1 - 1% BSA* in PBS, pH 7.2-7.4, 0.2 μm filtered.

IC Diluent #4** - 1 mM EDTA, 0.5% Triton X-100 in PBS, pH 7.2-7.4.

Note: Approximately 50 mL of this diluent is required to run the assay on one plate.

Lysis Buffer #12** - 1 mM EDTA, 0.005% Tween 20, 0.5% Triton X-100, 10 μ g/mL Leupeptin, 10 μ g/mL Pepstatin, 100 μ M PMSF, 3 μ g/mL Aprotinin in PBS, pH 7.2-7.4.

Substrate Solution - 1:1 mixture of Color Reagent A (H₂O₂) and Color Reagent B (Tetramethylbenzidine) (R&D Systems, Catalog # DY999).

Stop Solution - 2 N H₂SO₄ (R&D Systems, Catalog # DY994).

*The use of R&D Systems Reagent Diluent Concentrate 2 (Catalog # DY995) or Millipore Bovine Serum Albumin, Fraction V, Protease free (Catalog # 82-045) is recommended. All buffers containing BSA must be stored at 2-8° C.

**Sample Diluent Concentrate 1 (5X) (R&D Systems, Catalog # DYC001), prepared as described in the DYC001 insert.

Triton is a registered trademark of Union Carbide. Tween is a registered trademark of ICI Americas.

REAGENT PREPARATION

Bring all reagents to room temperature before use.

Total XIAP Capture Antibody (Part 840758) - Each vial contains 360 μ g/mL of mouse anti-human XIAP when reconstituted with 200 μ L of PBS. After reconstitution, store at 2-8° C for up to 30 days or aliquot and store at \leq -20° C in a manual defrost freezer or at \leq -70° C for up to 3 months.*

Total XIAP Detection Antibody (Part 840759) - Each vial contains 18 μ g/mL of biotinylated mouse anti-human XIAP when reconstituted with 1.0 mL of IC Diluent #1. After reconstitution, store at 2-8° C for up to 30 days or aliquot and store at \leq -20° C in a manual defrost freezer or at \leq -70° C for up to 3 months.*

Total XIAP Standard (Part 840760) - Each vial contains 80 ng/mL of recombinant human XIAP when reconstituted with 500 μ L of IC Diluent #4. Use a fresh standard for each assay. Use within one hour of reconstitution. A seven point standard curve using 2-fold serial dilutions in IC Diluent #4 and a high standard of 20,000 pg/mL is recommended.

Streptavidin-HRP (Part 890803) - Each vial contains 1.0 mL of streptavidin conjugated to horseradish-peroxidase. Store at 2-8° C. **DO NOT FREEZE.**

PREPARATION OF SAMPLES

Cell Lysates - Rinse cells two times with PBS, making sure to remove any remaining PBS after the second rinse. Solubilize cells at 1 x 10^7 cells/mL in Lysis Buffer #12 and allow samples to sit on ice for 15 minutes. Assay immediately or store at \leq -70° C. Before use, centrifuge samples at 2000 x g for 5 minutes and transfer the supernate to a clean test tube. Sample protein concentration may be quantified using a total protein assay. If needed, further dilutions should be made in IC Diluent #4.

^{*}Provided this is within the expiration date of the kit.

PRECAUTION

The Stop Solution suggested for use with this kit is an acidic solution. Wear protective gloves, clothing, eye, and face protection. Wash hands thoroughly after handling.

TECHNICAL HINTS AND LIMITATIONS

- This DuoSet IC ELISA should not be used beyond the expiration date on the kit label.
- Individual results may vary due to differences in technique, plasticware and water sources.
- It is important that the diluents selected for reconstitution and for dilution of the standard reflect the environment of the samples being measured. The diluent suggested in this protocol should be suitable for most cell lysates.
- The type of enzyme and substrate and the concentrations of capture/detection antibodies used can be varied to create an immunoassay with a different sensitivity and dynamic range. A basic understanding of immunoassay development is required for the successful use of these reagents in immunoassays.
- A thorough and consistent wash technique is essential for proper assay performance. Wash Buffer should be dispensed forcefully and removed completely from the wells by aspiration or decanting. Remove any remaining Wash Buffer by inverting the plate and blotting it against clean paper towels.
- Use a fresh reagent reservoir and pipette tips for each step.
- It is recommended that all standards and samples be assayed in duplicate.
- Avoid microbial contamination of reagents and buffers. This may interfere with the sensitivity of the assay. Buffers containing protein should be made under aseptic conditions and be stored at 2-8° C or prepared fresh daily.

GENERAL ELISA PROTOCOL

A plate layout is provided to record standards and samples assayed.

Plate Preparation

- 1. Dilute the Capture Antibody to a working concentration of 2.0 μ g/mL in PBS without carrier protein. Immediately coat a 96 well microplate with 100 μ L per well of the diluted Capture Antibody. Seal the plate and incubate overnight at room temperature.
- 2. Aspirate each well and wash with Wash Buffer, repeating the process two times for a total of 3 washes. Wash by filling each well with Wash Buffer (400 μ L) using a squirt bottle, manifold dispenser, or autowasher. Complete removal of liquid at each step is essential for good performance. After the last wash, remove any remaining Wash Buffer by aspirating or by inverting the plate and blotting it against clean paper towels.
- 3. Block plates by adding 300 μ L of Block Buffer to each well. Incubate at room temperature for 1-2 hours.
- 4. Repeat the aspiration/wash as in step 2. The plates are now ready for sample addition.

Assay Procedure

- 1. Add 100 μ L of sample or standards in IC Diluent #4 per well. Use IC Diluent #4 as the zero standard. Cover with a plate sealer and incubate 2 hours at room temperature.
 - **Note:** A seven point standard curve using 2-fold serial dilutions and a high standard of 20,000 pg/mL is recommended.
- 2. Repeat the aspiration/wash as in step 2 of Plate Preparation.
- 3. Dilute the Detection Antibody to a working concentration of 0.5 μ g/mL in IC Diluent #1 before use. Add 100 μ L of the diluted Detection Antibody to each well. Cover with a new plate sealer and incubate 2 hours at room temperature.
- 4. Repeat the aspiration/wash as in step 2 of Plate Preparation.
- 5. Immediately before use, dilute the Streptavidin-HRP to the working concentration specified on the vial label using IC Diluent #1. Add 100 μ L of the diluted Streptavidin-HRP to each well. Incubate for 20 minutes at room temperature. Avoid placing the plate in direct light.
- 6. Repeat the aspiration/wash as in step 2 of the Plate Preparation.
- 7. Add 100 μ L of Substrate Solution to each well. Incubate for 20 minutes at room temperature. Avoid placing the plate in direct light.
- 8. Add 50 μ L of Stop Solution to each well. Gently tap the plate to ensure thorough mixing.
- 9. Determine the optical density of each well immediately, using a microplate reader set to 450 nm. If wavelength correction is available, set to 540 nm or 570 nm. If wavelength correction is not available, subtract readings at 540 nm or 570 nm from the readings at 450 nm. This subtraction will correct for optical imperfections in the plate. Readings made directly at 450 nm without correction may be higher and less accurate.


CALCULATION OF RESULTS

Average the duplicate readings for each standard and sample then subtract the average zero standard optical density. Results may be normalized to total protein or cell number.

Create a standard curve by reducing the data using computer software capable of generating a four parameter logistic (4-PL) curve-fit. As an alternative, construct a standard curve by plotting the mean absorbance for each standard on the y-axis against the concentration on the x-axis and draw a best fit curve through the points on the graph. The data may be linearized by plotting the log of the XIAP concentrations versus the log of the O.D. and the best fit line can be determined by regression analysis. This procedure will produce an adequate but less precise fit of the data.

TYPICAL DATA

A standard curve should be generated for each set of samples assayed. The graph below represents typical data generated when using this Human Total XIAP DuoSet IC ELISA. The standard curve was calculated using a computer generated 4-PL curve-fit. This standard curve is for demonstration purposes only.

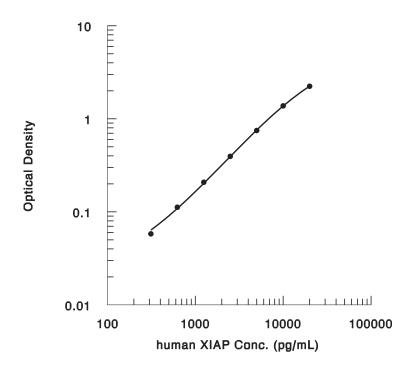

CALIBRATION

This DuoSet IC ELISA is calibrated against a highly purified *E. coli*-expressed recombinant human XIAP produced at R&D Systems.

SPECIFICITY

The Human Total XIAP Duoset IC ELISA is specific for XIAP. Specificity was demonstrated by Western blot analysis of protein bound by the capture antibody supplied in the kit.

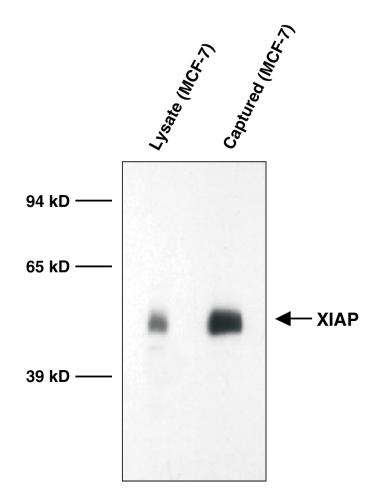


Figure 1: Lysates prepared from MCF-7 human breast cancer cells were incubated in wells coated with Total XIAP Capture Antibody. Unbound material was removed by washing and bound material was solubilized in SDS gel sample buffer. Captured proteins were electrophoresed, transferred to a PVDF membrane, and immunoblotted with Total XIAP DuoSet IC Detection Antibody. A single band corresponding to human XIAP was detected.

QUANTIFICATION

Amounts of human XIAP, as measured by this Duoset IC ELISA, are consistent with the relative amounts of human XIAP determined by qualitative Western blot analysis.


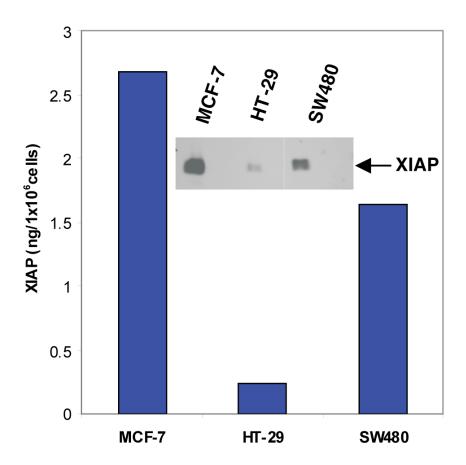


Figure 2: MCF-7 human breast cancer cells, HT-29 human colon adenocarcinoma cells, and SW480 human bladder cancer cells were lysed as described in the Preparation of Samples section. Human XIAP was measured with this DuoSet IC ELISA. The same lysates were immunoblotted (inset) with rabbit anti-human/mouse XIAP (R&D Systems, Catalog # AF822). The DuoSet IC ELISA results correlate well with the relative amounts of XIAP detected by Western blot.

PLATE LAYOUT

Use this plate layout as a record of standards and samples assayed.

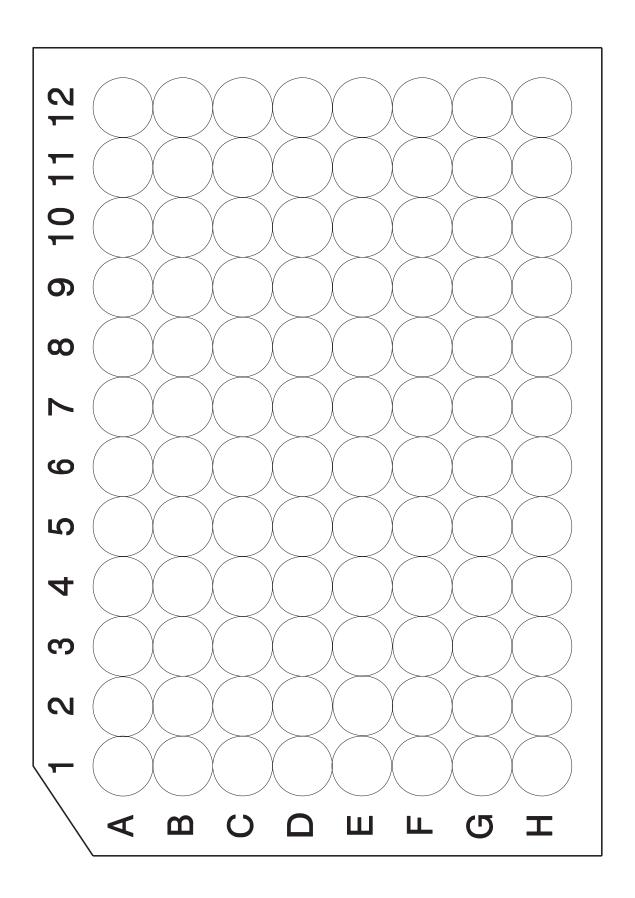

NOTES