

## DESCRIPTION


<b>Species Reactivity</b>	Human
<b>Specificity</b>	Detects human Endoglin/CD105.
<b>Source</b>	Monoclonal Mouse IgG <sub>1</sub> Clone # 166707
<b>Purification</b>	Protein A or G purified from hybridoma culture supernatant
<b>Immunogen</b>	Mouse myeloma cell line NS0-derived recombinant human Endoglin/CD105 Glu26-Gly586 Accession # Q5T9B9
<b>Conjugate</b>	Allophycocyanin Excitation Wavelength: 620-650 nm Emission Wavelength: 660-670 nm
<b>Formulation</b>	Supplied in a saline solution containing BSA and Sodium Azide. See Certificate of Analysis for details.  *Contains <0.1% Sodium Azide, which is not hazardous at this concentration according to GHS classifications. Refer to the Safety Data Sheet (SDS) for additional information and handling instructions.

## APPLICATIONS

**Please Note:** Optimal dilutions should be determined by each laboratory for each application. *General Protocols* are available in the *Technical Information* section on our website.

	Recommended Concentration	Sample
<b>Flow Cytometry</b>	10 $\mu$ L/10 <sup>6</sup> cells	See Below

## DATA

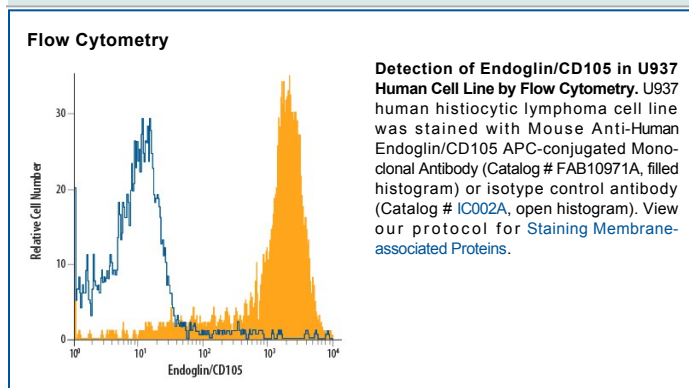

## PREPARATION AND STORAGE

<b>Shipping</b>	The product is shipped with polar packs. Upon receipt, store it immediately at the temperature recommended below.
<b>Stability &amp; Storage</b>	<b>Protect from light. Do not freeze.</b> <ul style="list-style-type: none"> <li>12 months from date of receipt, 2 to 8 °C as supplied.</li> </ul>

## BACKGROUND

Endoglin (CD105) is a 90 kDa type I transmembrane glycoprotein of the zona pellucida (ZP) family of proteins (1-3). Endoglin and betaglycan/TβR111 are type III receptors for TGF beta superfamily ligands, sharing 71% aa identity in the transmembrane (TM) and cytoplasmic domains. Endoglin is highly expressed on proliferating vascular endothelial cells, chondrocytes, and syncytiotrophoblasts of term placenta, with lower amounts on hematopoietic, mesenchymal and neural crest stem cells, activated monocytes, and lymphoid and myeloid leukemic cells (2-5). Human Endoglin cDNA encodes 658 amino acids (aa) including a 25 aa signal sequence, a 561 aa extracellular domain (ECD) with an orphan domain and a two-part ZP domain, a TM domain and a 47 aa cytoplasmic domain (1-3). An isoform with a 14 aa cytoplasmic domain (S-endoglin) can oppose effects of long (L) Endoglin (6, 7). The human Endoglin ECD shares 65-72% aa identity with mouse, rat, bovine, porcine and canine Endoglin. Endoglin homodimers interact with TGF-β1 and TGF-β3 (but not TGF-β2), but only after binding TβR111 (8). Similarly, they interact with activin-A and BMP-7 via activin type IIA or B receptors, and with BMP-2 via BMPR-1A/ALK-3 or BMPR-1B/ALK-6 (9). BMP-9, however, is reported to bind Endoglin directly (10). Endoglin modifies ligand-induced signaling in multiple ways. For example, expression of Endoglin can inhibit TGF-β1 signals but enhance BMP7 signals in the same myoblast cell line (11). In endothelial cells, Endoglin inhibits TβR1/ALK5, but enhances ALK1-mediated activation (12). Deletion of mouse Endoglin causes lethal vascular and cardiovascular defects, and human Endoglin haploinsufficiency can cause the vascular disorder, hereditary hemorrhagic telangiectasia type I (13, 14). These abnormalities confirm the essential function of Endoglin in differentiation of smooth muscle, angiogenesis, and neovascularization (2-4, 12-14). In preeclampsia of pregnancy, high levels of proteolytically generated soluble Endoglin and VEGF R1 (sFLT1), along with low placental growth factor (PlGF), are pathogenic due to antiangiogenic activity (15).

## References:

1. Gougos, A. and Letarte, M. (1990) *J. Biol. Chem.* **265**:8361.
2. ten Dijke, P. *et al.* (2008) *Angiogenesis* **11**:79.
3. Bernabeu, C. *et al.* (2007) *J. Cell. Biochem.* **102**:1375.
4. Mancini, M.L. *et al.* (2007) *Dev. Biol.* **308**:520.
5. Moody, J.L. *et al.* (2007) *Stem Cells* **25**:2809.
6. Velasco, S. *et al.* (2008) *J. Cell Sci.* **121**:913.
7. Perez-Gomez, E. *et al.* (2005) *Oncogene* **24**:4450.
8. Cheifetz, S. *et al.* (1992) *J. Biol. Chem.* **267**:19027.
9. Barbara, N.P. *et al.* (1999) *J. Biol. Chem.* **274**:584.
10. Scharpfenecker, M. *et al.* (2007) *J. Cell Sci.* **120**:964.
11. Scherner, O. *et al.* (2007) *J. Biol. Chem.* **282**:13934.
12. Pece-Barbara, N. *et al.* (2005) *J. Biol. Chem.* **280**:27800.
13. Arthur, H.M. *et al.* (2000) *Dev. Biol.* **217**:42.
14. Lebrin, F. and C.L. Mummery (2008) *Trends Cardiovasc. Med.* **18**:25.
15. Venkatesha, S. *et al.* (2006) *Nat. Med.* **12**:642.