

Magnetic Luminex® Performance Assay Human IL-15 High Sensitivity Kit

Catalog Number: LBHS247
Pack Size: 100 Tests

SPECIFICATIONS AND USE

Recommended Sample Types Microparticle Region Components

- Cell culture supernates, serum, EDTA plasma, and heparin plasma.
- Region-35
- Microparticle Concentrate (Part 894475) is supplied as a 50X concentrated stock (0.075 mL) with preservatives.
- Biotin-Antibody Concentrate (Part 894483) is supplied as a 100X concentrated stock solution (0.075 mL) with preservatives.

Other Supplies Required

 Magnetic Luminex Performance Assay Human High Sensitivity Cytokine Base Kit B (Catalog Number LBHS000).

Storage

- Store the unopened kit at 2-8 °C. Do not use past the expiration date on the label.
- Avoid freezing microparticles.
- Protect microparticles from light.


Instructions for Use

Refer to the Base Kit insert for the Magnetic Luminex Performance Assay procedure.

TYPICAL DATA

This human IL-15 standard curve is provided only for demonstration. A standard curve must be generated each time an assay is run, utilizing values from the Standard Value Card included in the Base Kit.

Note: When running cell culture supernate samples using Calibrator Diluent RD5K, a six-point standard curve (0.591-605 pg/mL) is recommended. When running serum/plasma samples using Calibrator Diluent RD6-65, a seven-point standard curve (0.591-2420 pg/mL) is recommended.

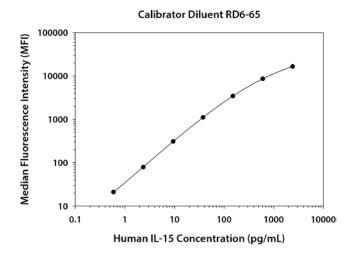

Standard	pg/mL	MFI	Average	Corrected	
Blank	0	23	24		
		24			
1	2420	16,375	16,471	16,447	
		16,566	10, 17 1	10,117	
2	605	8552	8622	0500	
Z	003	8691	8022	8598	
3	151	3438	3453	3429	
3	151	3467	3433		
4	37.8	1122	1132	1108	
		1141	1132		
5	9.45	331	333	309	
)		334	333		
6	2.36	101	103	79	
0		104	103		
7	0.591	44	45	21	
/		45	43	21	

752804.0

PRECISION

Intra-assay Precision (precision within an assay)

Three samples of known concentration were tested on one plate to assess precision within an assay.

Inter-assay Precision (precision between assays)

Three samples of known concentration were tested in separate assays to assess precision between assays.

Assays were performed by at least three technicians using two lots of components.

	Intra-assay Precision			Inter-assay Precision		
Sample	1	2	3	1	2	3
n	20	20	20	43	43	43
Mean (pg/mL)	4.92	59.0	617	4.79	57.7	541
Standard Deviation	0.127	0.812	8.25	0.607	6.39	68.7
% CV	2.6	1.4	1.3	12.7	11.1	12.7

RECOVERY & LINEARITY

Samples were spiked with human IL-15 and evaluated for recovery and were serially diluted to evaluate assay linearity.

Recovery					
Sample Type	Average % Recovery	Range (%)			
Cell culture supernates	81	74-101			
Serum	110	105-113			
EDTA plasma	99	92-109			
Heparin plasma	97	89-108			

Linearity						
		Cell culture supernates	Serum	EDTA Plasma	Heparin plasma	
1:2	Average % of Expected	110	95	95	99	
	Range (%)	108-112	93-98	94-99	96-102	
1:4	Average % of Expected	120	95	90	95	
	Range (%)	118-122	87-102	86-99	84-105	
1:8	Average % of Expected	126	92	86	95	
	Range (%)	125-127	78-98	78-103	79-110	

SENSITIVITY

All data were collected with assays run as a multiplex.

Twenty-one assays were evaluated, and the minimum detectable dose (MDD) of human IL-15 ranged from 0.034-0.167 pg/mL. The mean MDD was 0.065 pg/mL.

The MDD was determined by adding two standard deviations to the MFI of twenty zero standard replicates and calculating the corresponding concentration.

CORRELATION

This assay has been correlated to the Quantikine® ELISA Kit with a slope of 0.9-1.1 and an R² value greater than 0.9.

SPECIFICITY

Note: Refer to the base kit insert for a complete list of analytes tested for cross-reactivity and interference.

This assay recognizes natural and recombinant human IL-15.

TECHNICAL HINTS

- Protect the microparticles and streptavidin-PE from light at all times.
- Refer to the appropriate Base Kit Standard Value Card for reconstitution volume and values of the reconstituted standard.
- Diluted microparticles cannot be stored. Make a fresh dilution of microparticles each time the assay is run.
- The use of a magnetic device made to accommodate a microplate is necessary for washing.
- Discrepancies may exist in values obtained for the same analyte utilizing different technologies.

Magnetic Luminex Performance Assays afford the user the benefit of multianalyte analysis of cytokines in a complex sample. A single, multipurpose diluent for each sample type is used to optimize recovery, linearity, and reproducibility. Such a multipurpose, single diluent may not optimize any single analyte to the same degree that a unique diluent selected for analysis of that analyte can optimize conditions. Therefore, some performance characteristics may be more variable than those for assays designed specifically for single analyte analysis.

All trademarks and registered trademarks are the property of their respective owners.

752804.0 www.RnDSystems.com