



Magnetic Luminex® Performance Assay Human IL-22 High Sensitivity Kit

Catalog Number: LBHS5782

Pack Size: 100 Tests


SPECIFICATIONS AND USE

- Recommended Sample Types**
- Cell culture supernates, serum, EDTA plasma, and heparin plasma.
- Microparticle Region**
- Region-38
- Components**
- Microparticle Concentrate (Part 894478) is supplied as a 50X concentrated stock (0.075 mL) with preservatives.
 - Biotin-Antibody Concentrate (Part 894486) is supplied as a 100X concentrated stock solution (0.075 mL) with preservatives.
- Other Supplies Required**
- Magnetic Luminex Performance Assay Human High Sensitivity Cytokine Base Kit B (Catalog Number LBHS000).
- Storage**
- Store the unopened kit at 2-8 °C. Do not use past the expiration date on the label.
 - **Avoid freezing microparticles.**
 - **Protect microparticles from light.**
- Instructions for Use**
- Refer to the Base Kit insert for the Magnetic Luminex Performance Assay procedure.

TYPICAL DATA

This human IL-22 standard curve is provided only for demonstration. A standard curve must be generated each time an assay is run, utilizing values from the Standard Value Card included in the Base Kit.

Note: When running cell culture supernate samples using Calibrator Diluent RD5K, a six-point standard curve (1.64-1675 pg/mL) is recommended. When running serum/plasma samples using Calibrator Diluent RD6-65, a seven-point standard curve (1.64-6700 pg/mL) is recommended.

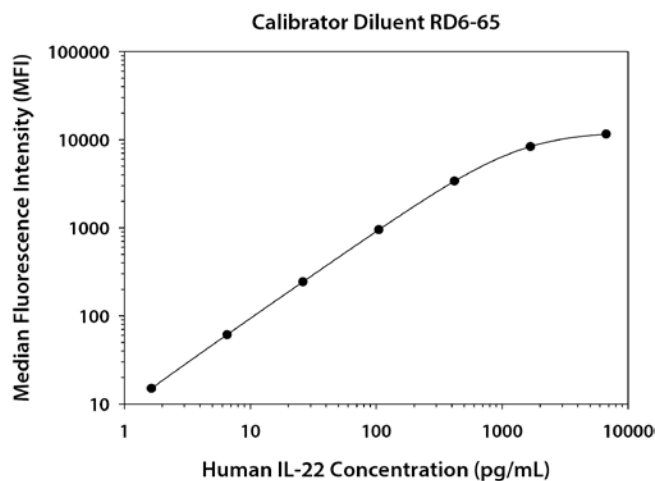

Standard	pg/mL	MFI	Average	Corrected
Blank	0	36 37	37	—
1	6700	11,544 11,666	11,605	11,568
2	1675	8361 8422	8392	8355
3	419	3405 3434	3419	3382
4	105	978 996	987	950
5	26.2	276 284	280	243
6	6.54	97 98	98	61
7	1.64	51 53	52	15

FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES.

R&D Systems, Inc.
1-800-343-7475

PRECISION

Intra-assay Precision (precision within an assay)

Three samples of known concentration were tested on one plate to assess precision within an assay.

Inter-assay Precision (precision between assays)

Three samples of known concentration were tested in separate assays to assess precision between assays.

Assays were performed by at least three technicians using two lots of components.

	Intra-assay Precision			Inter-assay Precision		
Sample	1	2	3	1	2	3
n	20	20	20	43	43	43
Mean (pg/mL)	6.93	72.2	663	7.76	84.5	709
Standard Deviation	0.209	1.13	11.0	1.4	11.2	104
% CV	3.0	1.6	1.7	18.0	13.3	14.7

RECOVERY & LINEARITY

Samples were spiked with human IL-22 and evaluated for recovery and were serially diluted to evaluate assay linearity.

Recovery			Linearity					
Sample Type	Average % Recovery	Range (%)			Cell culture supernates	Serum	EDTA Plasma	Heparin plasma
Cell culture supernates	100	93-108	1:2	Average % of Expected	100	107	101	107
				Range (%)	99-101	100-113	91-113	100-114
Serum	86	79-95	1:4	Average % of Expected	98	107	99	107
				Range (%)	97-98	95-121	83-116	95-120
EDTA plasma	84	70-100	1:8	Average % of Expected	102	103	94	104
				Range (%)	101-102	88-124	76-112	87-122
Heparin plasma	88	78-95						

SENSITIVITY

All data were collected with assays run as a multiplex.

Twenty-one assays were evaluated, and the minimum detectable dose (MDD) of human IL-22 ranged from 0.128-0.529 pg/mL. The mean MDD was 0.269 pg/mL.

The MDD was determined by adding two standard deviations to the MFI of twenty zero standard replicates and calculating the corresponding concentration.

CORRELATION

This assay has been correlated to the Quantikine® ELISA Kit with a slope of 0.9-1.1 and an R² value greater than 0.9.

SPECIFICITY

Note: Refer to the base kit insert for a complete list of analytes tested for cross-reactivity and interference.

This assay recognizes natural and recombinant human IL-22.

Recombinant human IL-1 RII interferes at concentrations > 25.0 ng/mL in this assay.

Recombinant human TNF RI interferes at concentrations > 25.0 ng/mL in this assay.

TECHNICAL HINTS

- Protect the microparticles and streptavidin-PE from light at all times.
- Refer to the appropriate Base Kit Standard Value Card for reconstitution volume and values of the reconstituted standard.
- Diluted microparticles cannot be stored. Make a fresh dilution of microparticles each time the assay is run.
- The use of a magnetic device made to accommodate a microplate is necessary for washing.
- Discrepancies may exist in values obtained for the same analyte utilizing different technologies.

Magnetic Luminex Performance Assays afford the user the benefit of multianalyte analysis of cytokines in a complex sample. A single, multipurpose diluent for each sample type is used to optimize recovery, linearity, and reproducibility. Such a multipurpose, single diluent may not optimize any single analyte to the same degree that a unique diluent selected for analysis of that analyte can optimize conditions. Therefore, some performance characteristics may be more variable than those for assays designed specifically for single analyte analysis.

All trademarks and registered trademarks are the property of their respective owners.