

Quantikine[®] ELISA

Mouse Dkk-1 Immunoassay

Catalog Number MKK100

For the quantitative determination of mouse Dickkopf related protein 1 (Dkk-1) concentrations in cell culture supernates, serum, and plasma.

Note: The standard reconstitution method has changed. Please read this package insert in its entirety before using this product.

This package insert must be read in its entirety before using this product.
For research use only. Not for use in diagnostic procedures.

TABLE OF CONTENTS

SECTION	PAGE
INTRODUCTION	1
PRINCIPLE OF THE ASSAY.....	2
LIMITATIONS OF THE PROCEDURE	2
TECHNICAL HINTS.....	2
MATERIALS PROVIDED & STORAGE CONDITIONS	3
OTHER SUPPLIES REQUIRED	3
PRECAUTIONS.....	4
SAMPLE COLLECTION & STORAGE	4
SAMPLE PREPARATION.....	4
REAGENT PREPARATION	5
ASSAY PROCEDURE	6
CALCULATION OF RESULTS.....	7
TYPICAL DATA.....	7
PRECISION	8
RECOVERY.....	8
LINEARITY.....	8
SENSITIVITY	9
CALIBRATION	9
SAMPLE VALUES.....	9
SPECIFICITY.....	9
REFERENCES.....	10

MANUFACTURED AND DISTRIBUTED BY:

USA & Canada | R&D Systems, Inc.

614 McKinley Place NE, Minneapolis, MN 55413, USA
TEL: (800) 343-7475 (612) 379-2956 FAX: (612) 656-4400
E-MAIL: info@RnDSystems.com

DISTRIBUTED BY:

UK & Europe | R&D Systems Europe, Ltd.

19 Barton Lane, Abingdon Science Park, Abingdon OX14 3NB, UK
TEL: +44 (0)1235 529449 FAX: +44 (0)1235 533420
E-MAIL: info@RnDSystems.co.uk

China | R&D Systems China Co., Ltd.

24A1 Hua Min Empire Plaza, 726 West Yan An Road, Shanghai PRC 200050
TEL: +86 (21) 52380373 FAX: +86 (21) 52371001
E-MAIL: info@RnDSystemsChina.com.cn

INTRODUCTION

Dickkopf related protein 1 (Dkk-1) is the founding member of the Dickkopf family of secreted proteins that includes Dkk-1, -2, -3, -4, and a related protein, Soggy (1, 2). Dickkopf, meaning "big head" or "stubborn", was discovered as an inducer of head formation in *Xenopus* (3). Dkk proteins contain two conserved cysteine-rich domains separated by a linker region (1-3). The C-terminal domain, which contains a colipase fold with a conserved pattern of ten cysteine residues, is necessary and sufficient for Wnt inhibition (4, 5). Mature mouse Dkk-1 is a 40 kDa glycosylated protein that shows 95%, 85%, 81%, and 83% amino acid (aa) sequence identity with rat, human, bovine, and rabbit Dkk-1, respectively. It also shows 42% and 38% aa sequence identity with mouse Dkk-2 and Dkk-4, respectively, with similarity mainly within the cysteine-rich domains.

Dkk-1 and Dkk-4 are well-documented antagonists of the canonical Wnt signaling pathway (1, 2). This pathway is activated by Wnt engagement of a receptor complex composed of the Frizzled proteins and one of two low-density lipoprotein receptor-related proteins, LRP-5 or LRP-6 (6). Dkk-1 antagonizes Wnt by forming ternary complexes of LRP-5/6 with Kremen-1 or Kremen-2 (7). Internalization of the Dkk-1/LRP-6/Kremen-2 complex downregulates Wnt signaling (6, 7). Dkk-1 has also been proposed to have Wnt-independent activity in some human cancer cell lines (8, 9). Dkk-1 is expressed throughout embryogenesis and antagonizes Wnt-7a during limb development, in developing neurons, keratinocytes, hair follicles, and the retina of the eye (10-14). Postnatal Dkk-1 is expressed mainly by osteoblasts and osteocytes (14).

The balance between Wnt signaling and Dkk-1 inhibition is critical for bone formation and homeostasis. Insufficient or excess Dkk-1 activity in bone results in increased or decreased bone density, respectively (14-16). High Dkk-1 expression has been shown and may be pathogenic in conditions where bone is eroded, such as rheumatoid arthritis, multiple myeloma, Paget's disease, and glucocorticoid-induced osteoporosis (17-22). Although the main phenotypes of experimental Dkk-1 deficiency are bone-related, it is important for regulating Wnt activity in other areas as well. Activity in the nervous system is indicated by the requirement of Dkk-1 expression for neural differentiation of mouse embryonic stem cells and for ischemic neuronal death (12, 23). Dkk-1 also regulates skin pigmentation and thickness by controlling Wnt signaling in melanocytes (13). Activation of Wnt by repression of Dkk-1 activity may be a factor in oncogenic transformation, for example, by the oncogene *c-myc* in mammary epithelial cell transformation or in human colon cancer (24, 25).

The Quantikine Mouse Dkk-1 immunoassay is a 4.5 hour solid-phase ELISA designed to measure Dkk-1 in cell culture supernates, serum, and plasma. It contains NS0-expressed recombinant mouse Dkk-1 and antibodies raised against the recombinant factor. This immunoassay has been shown to accurately quantitate the recombinant mouse Dkk-1. Results obtained using natural mouse Dkk-1 showed dose response curves that were parallel to the standard curves obtained using the Quantikine mouse kit standards. These results indicate that this kit can be used to determine relative mass values for natural mouse Dkk-1.

PRINCIPLE OF THE ASSAY

This assay employs the quantitative sandwich enzyme immunoassay technique. A monoclonal antibody specific for mouse Dkk-1 has been pre-coated onto a microplate. Standards, Control, and samples are pipetted into the wells and any mouse Dkk-1 present is bound by the immobilized antibody. After washing away any unbound substances, an enzyme-linked polyclonal antibody specific for mouse Dkk-1 is added to the wells. Following a wash to remove any unbound antibody-enzyme reagent, a substrate solution is added to the wells. The enzyme reaction yields a blue product that turns yellow when the Stop Solution is added. The intensity of the color measured is in proportion to the amount of mouse Dkk-1 bound in the initial step. The sample values are then read off the standard curve.

LIMITATIONS OF THE PROCEDURE

- FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES.
- The kit should not be used beyond the expiration date on the kit label.
- Do not mix or substitute reagents with those from other lots or sources.
- If samples generate values higher than the highest standard, further dilute the samples with Calibrator Diluent and repeat the assay.
- Any variation in standard diluent, operator, pipetting technique, washing technique, incubation time or temperature, and kit age can cause variation in binding.
- Variations in sample collection, processing, and storage may cause sample value differences.
- This assay is designed to eliminate interference by other factors present in biological samples. Until all factors have been tested in the Quantikine Immunoassay, the possibility of interference cannot be excluded.

TECHNICAL HINTS

- When mixing or reconstituting protein solutions, always avoid foaming.
- To avoid cross-contamination, change pipette tips between additions of each standard level, between sample additions, and between reagent additions. Also, use separate reservoirs for each reagent.
- For best results, pipette reagents and samples into the center of each well.
- The samples must be added to the plate within 15 minutes.
- To ensure accurate results, proper adhesion of plate sealers during incubation steps is necessary.
- Substrate Solution should remain colorless until added to the plate. Keep Substrate Solution protected from light. Substrate Solution should change from colorless to gradations of blue.
- Stop Solution should be added to the plate in the same order as the Substrate Solution. The color developed in the wells will turn from blue to yellow upon addition of the Stop Solution.

MATERIALS PROVIDED & STORAGE CONDITIONS

Store the unopened kit at 2-8 °C. Do not use past kit expiration date.

PART	PART #	DESCRIPTION	STORAGE OF OPENED/ RECONSTITUTED MATERIAL
Mouse Dkk-1 Microplate	893325	96 well polystyrene microplate (12 strips of 8 wells) coated with a monoclonal antibody specific for mouse Dkk-1.	Return unused wells to the foil pouch containing the desiccant pack. Reseal along entire edge of the zip-seal. May be stored for up to 1 month at 2-8 °C.*
Mouse Dkk-1 Standard	893327	2 vials of recombinant mouse Dkk-1 in a buffered protein base with preservatives; lyophilized. <i>Refer to the vial label for reconstitution volume.</i>	Discard after use. Use a fresh Standard and Control for each assay.
Mouse Dkk-1 Control	893328	2 vials of recombinant mouse Dkk-1 in a buffered protein base with preservatives; lyophilized. The assay value of the Control should be within the range specified on the label.	
Mouse Dkk-1 Conjugate	893326	12 mL of a polyclonal antibody specific for mouse Dkk-1 conjugated to horseradish peroxidase with preservatives.	May be stored for up to 1 month at 2-8 °C.*
Assay Diluent RD1-21	895215	12 mL of a buffered protein base with preservatives.	
Calibrator Diluent RD5-26 Concentrate	895525	21 mL of a concentrated buffered protein base with preservatives.	
Wash Buffer Concentrate	895003	21 mL of a 25-fold concentrated solution of buffered surfactant with preservative. <i>May turn yellow over time.</i>	
Color Reagent A	895000	12 mL of stabilized hydrogen peroxide.	
Color Reagent B	895001	12 mL of stabilized chromogen (tetramethylbenzidine).	
Stop Solution	895174	23 mL of diluted hydrochloric acid.	
Plate Sealers	N/A	4 adhesive strips.	

* Provided this is within the expiration date of the kit.

OTHER SUPPLIES REQUIRED

- Microplate reader capable of measuring absorbance at 450 nm, with the correction wavelength set at 540 nm or 570 nm.
- Pipettes and pipette tips.
- Deionized or distilled water.
- Squirt bottle, manifold dispenser, or automated microplate washer.
- 100 mL and 500 mL graduated cylinders.
- **Polypropylene** test tubes for dilution of standards and samples.

PRECAUTIONS

The Stop Solution provided with this kit is an acid solution.

Some components in this kit contain ProClin® which may cause an allergic skin reaction. Avoid breathing mist.

Color Reagent B may cause skin, eye, and respiratory irritation. Avoid breathing fumes.

Wear protective gloves, clothing, eye, and face protection. Wash hands thoroughly after handling. Please refer to the MSDS on our website prior to use.

SAMPLE COLLECTION & STORAGE

The sample collection and storage conditions listed below are intended as general guidelines. Sample stability has not been evaluated.

Cell Culture Supernates - Remove particulates by centrifugation. Assay immediately or aliquot and store samples at ≤ -20 °C. Avoid repeated freeze-thaw cycles.

Serum - Allow blood samples to clot for 2 hours at room temperature before centrifuging for 20 minutes at 2000 x g. Remove serum and assay immediately or aliquot and store samples at ≤ -20 °C. Avoid repeated freeze-thaw cycles.

Plasma - Collect plasma using EDTA or heparin as an anticoagulant. Centrifuge for 20 minutes at 2000 x g within 30 minutes of collection. Assay immediately or aliquot and store samples at ≤ -20 °C. Avoid repeated freeze-thaw cycles.

Note: *Citrate plasma has not been validated for use in this assay.*

SAMPLE PREPARATION

Serum and plasma samples require a 10-fold dilution prior to assay. A suggested 10-fold dilution is 20 μ L of sample + 180 μ L of Calibrator Diluent RD5-26 (diluted 1:4).*

* See Reagent Preparation section.

All trademarks and registered trademarks are the property of their respective owners.

REAGENT PREPARATION

Bring all reagents to room temperature before use.

Mouse Dkk-1 Control - Reconstitute the Control with 1.0 mL of deionized or distilled water. Mix thoroughly. Assay the Control undiluted.


Wash Buffer - If crystals have formed in the concentrate, warm to room temperature and mix gently until the crystals have completely dissolved. Add 20 mL of Wash Buffer Concentrate to deionized or distilled water to prepare 500 mL of Wash Buffer.

Calibrator Diluent RD5-26 (diluted 1:4) - Add 20 mL of Calibrator Diluent RD5-26 Concentrate to 60 mL of deionized or distilled water to prepare 80 mL of Calibrator Diluent RD5-26 (diluted 1:4).

Substrate Solution - Color Reagents A and B should be mixed together in equal volumes within 15 minutes of use. Protect from light. 100 μ L of the resultant mixture is required per well.

Mouse Dkk-1 Standard - Refer to the vial label for reconstitution volume. Reconstitute the Mouse Dkk-1 Standard with Calibrator Diluent RD5-26 (diluted 1:4). Do not substitute other diluents. This reconstitution produces a stock solution of 4000 pg/mL. Allow the stock solution to sit for a minimum of 5 minutes with gentle mixing prior to making dilutions.

Use polypropylene tubes. Pipette 200 μ L of Calibrator Diluent RD5-26 (diluted 1:4) into each tube. Use the stock solution to produce a 2-fold dilution series (below). Mix each tube gently but thoroughly before the next transfer. The undiluted Mouse Dkk-1 Standard (4000 pg/mL) serves as the high standard. Calibrator Diluent RD5-26 (diluted 1:4) serves as the zero standard (0 pg/mL).

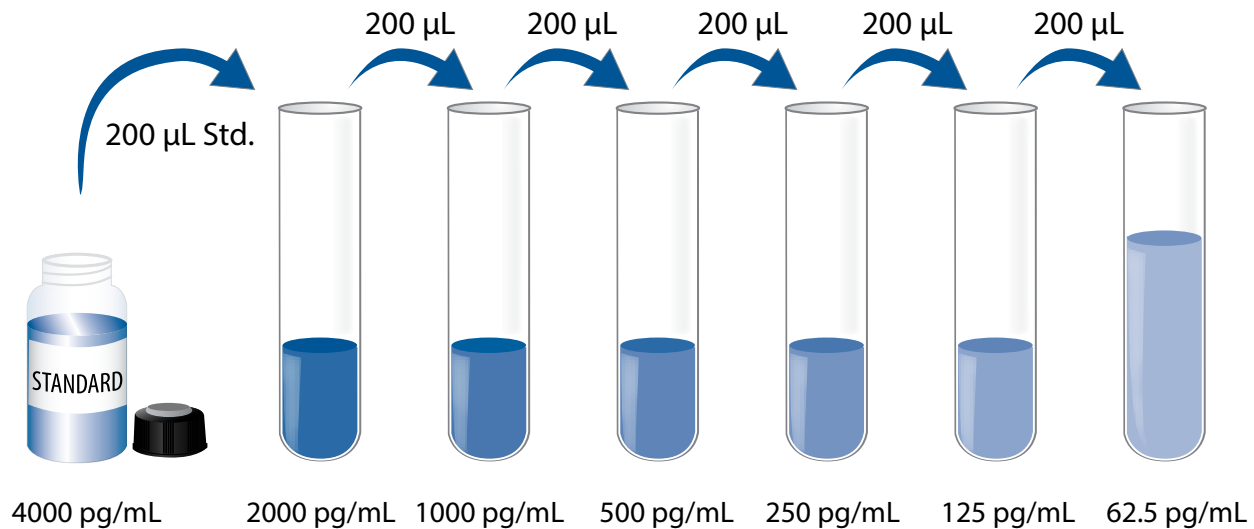

ASSAY PROCEDURE

Bring all reagents and samples to room temperature before use. It is recommended that all standards, Control, and samples be assayed in duplicate.

1. Prepare all reagents, standard dilutions, Control, and samples as directed in the previous sections.
2. Remove excess microplate strips from the plate frame, return them to the foil pouch containing the desiccant pack, and reseal.
3. Add 50 μL of Assay Diluent RD1-21 to each well.
4. Add 50 μL of Standard, Control, or sample* per well. **Samples must be added to the plate within 15 minutes.** Cover with the adhesive strip provided. Incubate for 2 hours at room temperature.
5. Aspirate each well and wash, repeating the process four times for a total of five washes. Wash by filling each well with Wash Buffer (400 μL) using a squirt bottle, manifold dispenser, or autowasher. Complete removal of liquid at each step is essential to good performance. After the last wash, remove any remaining Wash Buffer by aspirating or decanting. Invert the plate and blot it against clean paper towels.
6. Add 100 μL of Mouse Dkk-1 Conjugate to each well. Cover with a new adhesive strip. Incubate for 2 hours at room temperature.
7. Repeat the aspiration/wash as in step 5.
8. Add 100 μL of Substrate Solution to each well. Incubate for 30 minutes at room temperature. **Protect from light.**
9. Add 100 μL of Stop Solution to each well. Gently tap the plate to ensure thorough mixing.
10. Determine the optical density of each well within 30 minutes, using a microplate reader set to 450 nm. If wavelength correction is available, set to 540 nm or 570 nm. If wavelength correction is not available, subtract readings at 540 nm or 570 nm from the readings at 450 nm. This subtraction will correct for optical imperfections in the plate. Readings made directly at 450 nm without correction may be higher and less accurate.

*Samples may require dilution. See the Sample Preparation section.

CALCULATION OF RESULTS


Average the duplicate readings for each standard, Control, and sample and subtract the average zero standard optical density (O.D.).

Create a standard curve by reducing the data using computer software capable of generating a four parameter logistic (4-PL) curve-fit. As an alternative, construct a standard curve by plotting the mean absorbance for each standard on the y-axis against the concentration on the x-axis and draw a best fit curve through the points on the graph. The data may be linearized by plotting the log of the mouse Dkk-1 concentrations versus the log of the O.D. and the best fit line can be determined by regression analysis. This procedure will produce an adequate but less precise fit of the data.

If samples have been diluted, the concentration read from the standard curve must be multiplied by the dilution factor.

TYPICAL DATA

This standard curve is provided for demonstration only. A standard curve should be generated for each set of samples assayed.

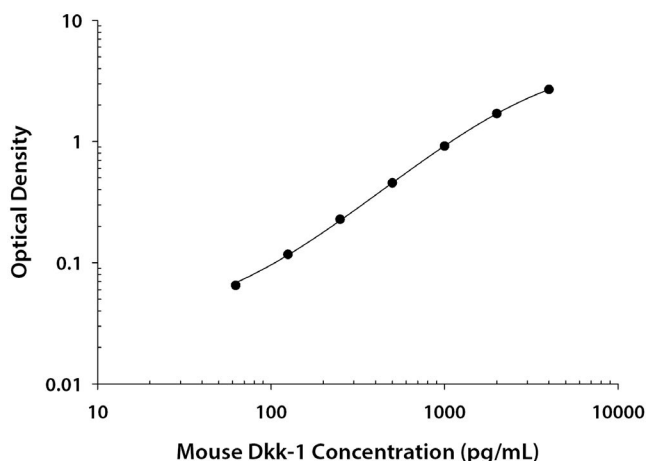

(pg/mL)	O.D.	Average	Corrected
0	0.070 0.071	0.071	—
62.5	0.135 0.136	0.136	0.065
125	0.186 0.189	0.188	0.117
250	0.297 0.300	0.299	0.228
500	0.524 0.528	0.526	0.455
1000	0.989 0.989	0.989	0.918
2000	1.750 1.789	1.770	1.699
4000	2.723 2.791	2.757	2.686

PRECISION

Intra-assay Precision (Precision within an assay)

Three samples of known concentration were tested twenty times on one plate to assess intra-assay precision.

Inter-assay Precision (Precision between assays)

Three samples of known concentration were tested in forty separate assays to assess inter-assay precision. Assays were performed by at least three technicians using two lots of kit components.

Sample	Intra-Assay Precision			Inter-Assay Precision		
	1	2	3	1	2	3
n	20	20	20	40	40	40
Mean (pg/mL)	230	478	1487	219	481	1492
Standard deviation	20.1	29.7	80.8	20.8	44.7	144
CV (%)	8.7	6.2	5.4	9.5	9.3	9.7

RECOVERY

The recovery of mouse Dkk-1 levels throughout the range of the assay in various matrices was evaluated.

Sample Type	Average % Recovery	Range
Cell culture supernates (n=4)	99	88-107%
Serum* (n=4)	113	100-120%
EDTA plasma* (n=4)	109	93-117%
Heparin plasma* (n=4)	113	104-119%

*Samples were diluted prior to assay.

LINEARITY

To assess the linearity of the assay, samples containing and/or spiked with high concentrations of mouse Dkk-1 were diluted with Calibrator Diluent to produce samples with values within the dynamic range of the assay.

		Cell culture supernates (n=4)	Serum* (n=4)	EDTA plasma* (n=4)	Heparin plasma* (n=4)
1:2	Average % of Expected	95	97	101	101
	Range (%)	93-99	93-104	95-107	95-109
1:4	Average % of Expected	91	96	101	101
	Range (%)	85-95	89-110	90-110	89-111
1:8	Average % of Expected	89	102	105	98
	Range (%)	81-97	94-116	93-114	87-108
1:16	Average % of Expected	88	98	100	—
	Range (%)	81-97	84-117	90-120	—

*Samples were diluted prior to assay as directed in the Sample Preparation section.

SENSITIVITY

Seventy-nine assays were evaluated and the minimum detectable dose (MDD) of mouse Dkk-1 ranged from 5.1-29.6 pg/mL. The mean MDD was 13.4 pg/mL.

The MDD was determined by adding two standard deviations to the mean optical density value of twenty zero standard replicates and calculating the corresponding concentration.

CALIBRATION

This immunoassay is calibrated against a highly purified NS0-expressed recombinant mouse Dkk-1 produced at R&D Systems.

SAMPLE VALUES

Serum/Plasma - Samples were evaluated for the presence of mouse Dkk-1 in this assay.

Sample Type	Mean (pg/mL)	Range (pg/mL)	Standard Deviation (pg/mL)
Serum (n=19)	13,893	1610-27,564	7579
EDTA plasma (n=19)	10,069	1616-24,500	7394
Heparin plasma (n=19)	4522	1446-9458	2549

SPECIFICITY

This assay recognizes natural and recombinant mouse Dkk-1.

The factors listed below were prepared at 40 ng/mL in Calibrator Diluent and assayed for cross-reactivity. Preparations of the following factors at 40 ng/mL in a mid-range mouse Dkk-1 control were assayed for interference. No significant cross-reactivity or interference was observed.

Recombinant mouse:

Dkk-2
Dkk-3
Dkk-4
Dkk-4 (R67S/R70S mutant)
Kremen-1
Kremen-2
LRP-6

This assay has not been validated for rat samples, however, recombinant rat Dkk-1 cross-reacts approximately 94%.

Recombinant human Dkk-1 cross-reacts approximately 26% in this assay.

REFERENCES

1. Krupnik, V.E. *et al.* (1999) *Gene* **238**:301.
2. Niehrs, C. (2006) *Oncogene* **25**:7469.
3. Glinka, A. *et al.* (1998) *Nature* **389**:517.
4. Bullock, C.M. *et al.* (2004) *Mol. Pharmacol.* **65**:582.
5. Brott, B.K. and S.Y. Sokol (2002) *Mol. Cell. Biol.* **22**:6100.
6. Mao, B. *et al.* (2001) *Nature* **411**:321.
7. Mao, B. *et al.* (2002) *Nature* **417**:664.
8. Lee, A.Y. *et al.* (2004) *Biochem. Biophys. Res. Commun.* **323**:1246.
9. Peng, S. *et al.* (2006) *Biochem. Biophys. Res. Commun.* **350**:641.
10. Kemp, C. *et al.* (2005) *Dev. Dyn.* **233**:1064.
11. Adamska, M. *et al.* (2004) *Dev. Biol.* **272**:134.
12. Verani, R. *et al.* (2007) *J. Neurochem.* **100**:242.
13. Yamaguchi, Y. *et al.* (2007) *J. Invest. Dermatol.* **127**:1217.
14. Li, J. *et al.* (2006) *Bone* **36**:754.
15. Morvan, F. *et al.* (2006) *J. Bone Miner. Res.* **21**:934.
16. MacDonald, B.T. *et al.* (2007) *Bone* **41**:331.
17. Diarra, D. *et al.* (2007) *Nature Med.* **13**:156.
18. Tian, E. *et al.* (2003) *N. Engl. J. Med.* **349**:2483.
19. Naot, D. *et al.* (2007) *J. Bone Miner. Res.* **22**:298.
20. Ohnaka, K. *et al.* (2005) *Biochem. Biophys. Res. Commun.* **329**:177.
21. Yaccoby, S. *et al.* (2007) *Blood* **109**:2106.
22. Giuliani, N. *et al.* (2007) *Can. Res.* **67**:7665.
23. Cappuccio, I. *et al.* (2005) *J. Neurosci.* **25**:2647.
24. Cowling, V.H. *et al.* (2007) *Mol. Cell. Biol.* **27**:5135.
25. Gonzales-Sancho, J.M. *et al.* (2005) *Oncogene* **24**:1098.